
KULLANIM KILAVUZU

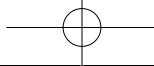TÜRKÇE VER.

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

ENGLISH VER.

This product is the beauty device for lip care with LED.
Please read this user manual carefully before operation and then use the product properly and safely.

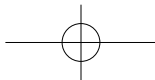Product Contents

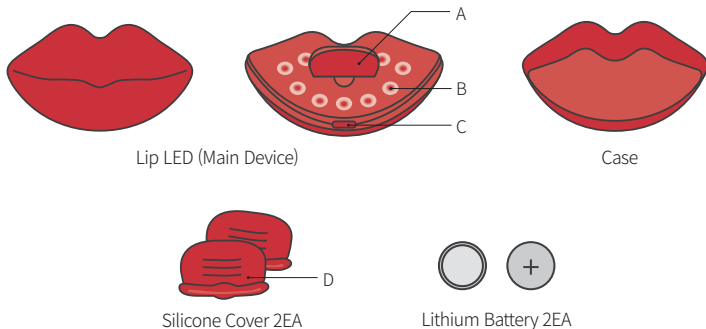A	Fixed Bar	The area is bitten with the teeth to hold
B	LED	RED 623nm(normal)

C	Power Switch	Turn on/off
D	Silicone Cover	Cover the fixed bar to use

* The contents may be different from the pictures by the model and the company regulations.

How to Use

- ① Apply the exclusive solution for BIBILIP LED (HONEY COCO TOTAL CARE) that comes with the device on your lip. (keep the solution on your lips for 10mins)


-After cleansing and using skin care products, apply enough amount of the solution. (keep the solution on your lips for 10mins)

* If you want to use your own lip moisturizing solutions that you already have, please apply enough amount of the solution after cleansing and using skin care products. (Keep the solution for 10mins)


* To use the Honey Coco Total Care solution only without the device, apply enough amount of the solution before going to sleep.

- ② After removing the case, insert the lithium battery that is provided into the device. (Refer to the page 5 of how to install the battery)
- ③ Cover the Silicone cover on the fixed bar

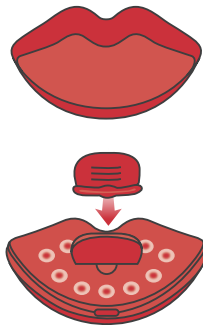
2/3

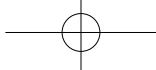- ④ Press the power switch to turn on the device.

- ⑤ Bite the fixed bar that covered with silicone cover slightly and then adhere the front of the lips to the back side of the device gently.

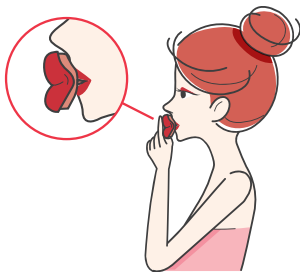
- ⑥ The device will be turn off automatically after 3mins. (To turn off the device while operating, press power switch)

- ⑦ After finishing the use of the device, separate the silicone cover from the device, wash the cover in the running water, and then wipe the back side of the device with wet-tissue.

*** BIBILIP Honey Coco Dual Moisturizing
(Use frequently)**



1. After cleansing and using skin care products, apply enough amount of the solution.
2. Apply your favorite color of lipstick
3. Use the solution frequently when you feel dry.

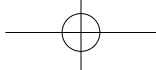How to Replace Lithium Battery


Flashing LED means low-battery and LED is not appeared even if press the power switch when the battery is depleted.

- ① Separate the fixed bar from the device by revolving the fixed bar 90degree clockwise.
- ② Take apart the cover.
- ③ Insert two 3V Lithium battery (CR2032) with correct polarity. (Check up on +,- direction)
- ④ Recombine the cover to the device with correct direction.
- ⑤ Put the fixed bar in the LED hole and turn the bar 90-degree counterclockwise to make them together.

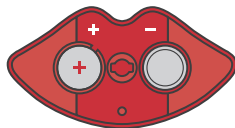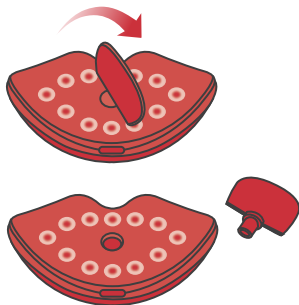*The battery can be operated more than 100 times based on using 3 minutes once day.

4/5

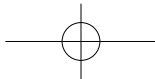
Product Specifications

Product Name	BIBILIP
Product No	BBL-100
Electric Power	DC 6V (3V Lithium Battery (CR2032) 2EA)
Dimension	12x9x8.5cm / 160g (Full package)
LED Wavelength	RED 623nm(normal)
Manufacturer	GTG Wellness
Address	A-902, Bundang Suji U-Tower, 767(Dongcheon-dong), Sinsu-ro, Suji-gu, Yongin-si, Gyeonggi-do, Korea, 16827
A/S Contact No	031-1600-2614
Manufactured Country	South Korea

Precaution

Warning

01. Read carefully the operation instruction and the precaution before operating the device.
02. The device is prohibited to use when you are exposed to sunlight and have sun-burn on your skin.
03. The recommendation of the device is once in daily. If you use the device extra, it may occur skin troubles.
04. If you have the symptoms as below, please use the device after the professional medical counsel.
 - Having skin trouble or medical treatment (allergic skin, dermatitis)
 - Being treated with dental care (implant, dentures, gum treatment)

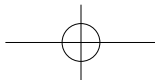
- Having epilepsy symptoms
 - Getting infectious diseases
 - Being oversensitive to light
 - Having sensitive skin, abnormal conditions, other diseases, or physical disorder
05. If you have skin trouble (skin redness, bleeding, itchy, etc.) during the use of the device, please stop use the device immediately and have medical counselling.
 06. Lip LED device can occur a slight temperatures while using the device. (If high temperature is occurred, turn off the device and then call the A/S center.)
 07. The device is not recommended to use for pregnant, menstrual, breast-feeding women.
 08. The device is not recommended for children under the 16-year-old.
 09. Use the device during driving, walking, or others is not allowed.
 10. When you feel dizzy or physical disorder during use of the device, please use the device after medical counselling.
 11. Staring the LED light directly may cause eyesight loss. Please avoid staring the LED light.

Precautions of Lithium Battery

01. Insulate the battery with the tape for storage of the battery.
02. Do not store the Lithium battery with other battery or metal goods.
03. For the battery disposal, insulate the battery with the tape and inquiry the district office about the battery disposal. If you do not insulate the battery properly, the fire by heating will be possible to occur.
04. As Lithium battery is consumable goods, the battery life is reduced by using the device.
05. Store the battery in the place where the children can not touch.
06. Swallowing the Lithium battery may cause serious injury or death. Take the treatment from doctors immediately.
07. The explosion or Leakage is possible by heating, disassembling, compressing, recharging, exposing the fire, or inserting into the product in the wrong way with the Lithium battery.

6/7

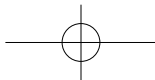Precautions

01. As this product is not waterproof, be careful that the product is not soaked in the water.
02. Operating the product with wet-hand may be the reason of malfunction.
03. Be careful not to damage the product. Damaged product can cause malfunction of the product.
04. Do not use the damaged product.
05. Do not disassemble/remodel the product on your own.
06. Be careful that using the device in dry conditions may cause static electricity.
07. Avoid to operate or store the product near fire or flammable materials.
08. Do not use the product for other purposes except the purposes described in this user manual.

Storage

01. Avoid to keep the device in the place with direct sunlight, high temperature or humidity.
02. After use of the device, separate the silicone cover from the device and wash the cover with the running water and then let them dry.
03. After use of the device, wipe the inside(back side) of the device with wet-tissue, wipe the device with dry clothes and then put it in the case for the storage. Do not use the strong chemical goods(ex. Alcohol) that can damage the surface of the device.
04. Due to the possibility of that the device can be damaged by external collision or pressure, keep the device in the safe place.
05. Avoid to keep the device in the place where children or pets can touch.
06. Do not put the device near volatile material.

Product Warranty

- This product is passed the inspection by our quality control standard. If the product is out of order or malfunction, the product will be repaired as following this warranty certificate. (This warranty certificate is not reissued.)
- Apply the A/S service to the stores where the device is purchased or the main office with the warranty certificate.

Product Name [Model No.]		Warranty period
Purchase Date		
Manufacture No.		
Customer		
Customer Address		
Customer Contact No.		

Warranty Service

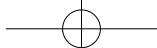
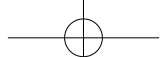
01. The period of warranty service is 1 year from the purchase date.
02. If there is manufacturer's defect or product malfunction in the warranty period, the product will be compensated under the free repair or consumer damage compensation regulations.
03. Under the following the conditions, the product is repaired with extra payment unless during the warranty period.
 - Failure due to the users' misuse or careless
 - Malfunction or Damages due to the product repair by others that are not selected A/S Center
 - Use the consumable goods, optional goods that are not provided with the device (Inquiry to the company for selected consumable goods and optional goods)
 - Natural Calamities (Breakdown by fire, saltwater, flood, etc.)
 - Users' Intentional purpose of damage or negligence without following the manual instructions

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