Turkish Ver.
English Ver.

LEBODY

FACE

USER MANUAL

LEBODY

— FACE —

English Ver.

LEBODY

— FACE —

WHAT IS LEBODY FACE?

Lebody Face is a home care device that gives vital, elastic and lifting effects to skin cells using micro current. Lebody Face makes amazing skin changes with only short time of use.

CONTENTS


- ① OVERVIEW**
PACKAGE COMPONENTS
- ② PRODUCT DESCRIPTION**
LEBODY FACE
- ③ HOW TO USE**
ORDER OF USING LEBODY FACE
HOW TO USE LEBODY FACE
LEBODY FACE LED DISPLAY
BATTERY CHARGING STATUS
- ④ PRECAUTIONS**
- ⑤ WARNINGS / PRODUCT CARE AFTER USE**
- ⑥ PRODUCT SPECIFICATIONS**
- ⑦ WARRANTY**

1 OVERVIEW

• PACKAGE COMPONENTS

The LEBODY FACE package consists of the following items

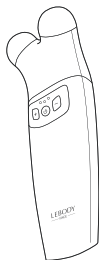LEBODY FACE

Charging Cradle

USB Charger Cable

User Manual


Pouch

※ The items may not look exactly as shown in the image, depending on the model and company regulations.

2 PRODUCT DESCRIPTION

• LEBODY FACE

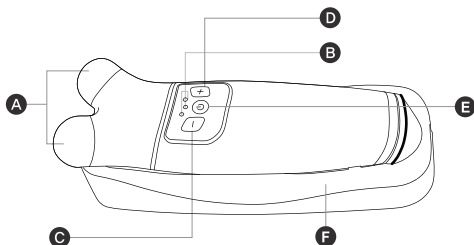No	NAME	FUNCTION
A	Microcurrent output unit (Faceball)	The area where the microcurrent is generated
B	LED display	The area where the mode and intensity are indicated with the LED lamp ● I mode ● L mode
C	Microcurrent intensity controller: down(-) / mode switch	A button for decreasing the output intensity. (Press and hold the button for 2~3 sec to switch between the I and L modes)
D	Microcurrent intensity controller: up(+) / mode switch	A button for increasing the output intensity. (Press and hold the button for 2~3 sec to switch between the I and L modes)
E	Power on/off	Press and hold for 2~3 sec to turn the power on/off
F	Charging Cradle	A holder for the LEBODY FACE during the charging process

● **I mode**

This mode promotes the absorption of facial creams. Apply Lebody renewal up cream, and use the device to massage the entire face for 5 to 10 min.

● **L mode**

This mode enhances skin elasticity with microcurrent stimulation. Apply Renewal face water gel, and use the device to massage the entire face for 5 to 10 min.

User Guidelines for LEBODY FACE

③ HOW TO USE

① Cleanse the face

Cleanse your face and remove water by towel. For men, shaving should be made before using the device.

② Apply the product

Take Renewal Up Cream or a Face Water Gel depending on your purpose of use and apply it evenly all over the face. (If you want to use the **I mode** and **L mode** back to back, Apply a Renewal Up Cream first at **I mode** and apply Face Water Gel at **L mode**.)

Sold Separately

● **Renewal up cream for I mode**

Renewal up cream for **I mode** (Briefly Up cream hereinafter)


● **Renewal Face water gel for L mode**

Renewal face water gel for **L mode** (Briefly Water gel hereinafter)

③ Turn on the power ④ Choose a mode ⑤ Adjust the intensity

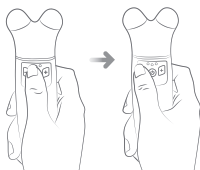Press and hold the power button for 2~3 sec to turn on the device.

Press and hold the [-/+] button for 2~3 sec to choose the mode you wish to use.

Press the [-/+] button to adjust the intensity of the microcurrent.

● **I mode** (Default mode)

● **L mode**

Recommended duration for each mode: 5~10 min

⑥ Use the device

Be sure to contact both Faceballs together on the face.

LEBODY FACE will be automatically turned off 10 minutes.

User Guidelines for the device

Proceed in the order and direction shown in the illustrations. Glide the device upward in a lifting motion. (Make sure the Faceballs are kept in close contact with the skin.) If you have problem-prone skin, make sure to cleanse your face again with warm water after using the device.

Starting point


Switch to a different area when you hear the beep.

4 seconds

Move the device from the starting point toward the direction of the arrow.

STEP 1

Glide the device as following illustration. Relocate the Faceballs to the next area every four seconds when you hear the beep. Be sure to contact both Faceballs together on the face.

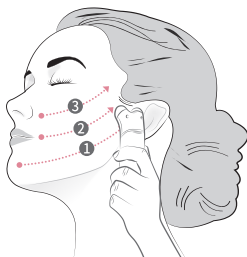STEP 2

※ Proceed in the same manner for the opposite side.

• LEBODY FACE LED DISPLAY

- **I mode:** Blue LED lamp is on
(default mode when the power is turned on)
- **L mode:** Red LED lamp is on
- Low battery: Red LED flickers five times and the power turns off
- Charging: Blue LED lamp turns on from left to right
- Fully charged: Blue LED lamp flickers

Low battery	Charging	Fully charged
Red LED flickers five times before the power turns off	Blue LED lamp turns on from left to right	Blue LED lamp flickers

• BATTERY CHARGING STATUS

- Turn on LEBODY FACE and place it on the charging cradle. The blue LED lamp will flicker from left to right to indicate that the device is being charged. When you place the device on the charging cradle, and the blue LED flickers, it means that the battery is fully charged.
- The device is released from the factory with a 30% charged battery. Please charge the device before using it for the first time. If the battery is significantly low due to natural discharge, it may take quite time to fully charge it.
- It takes around 2 and half hours to fully charge the battery. The battery will last for at least 7 hours at I mode and 5 and half hours at L mode, on average, when the device is used at level 3.
- If the device will not be used for a long time, make sure to fully charge the device before storing it away and fully charge it once a month as a way to prevent the battery from becoming overdischarged.
- If the device is not recharged for more than a month, it may not turn on because the battery has been overdischarged. In this case, use the device after recharging the battery.
- If the power does not turn on, if the power turns off immediately after the device has been turned on, or if the red LED lamp flickers, it means that the battery has been discharged. In this case, recharge the battery before use. If the battery has been completely discharged, it may take slightly longer than normal to fully charge it.

4 PRECAUTIONS

PRECAUTIONS

- Make sure to read and familiarize yourself with the instructions provided in the User Manual before using the device.
- This product is for facial care only. We cannot be held liable for any harm or damage resulting from the use of the device on any other areas of the body, aside from the face, or from inappropriate usage.
- Make sure to take off any jewelries (e.g. necklace, earrings, rings) prior to use.
- Do not use LEBODY FACE if you are allergic to metal or have sensitive skin.
- Do not store or use the device in humid areas, or in areas exposed to direct sunlight, or near a heating apparatus.
- When handling the device, do not drop it or subject it to strong impact.
- Make sure to store the device in an area that is out of reach of children.
- This product has not been waterproofed. It should not be used in water or washed with water but the external parts may be cleaned with wet wipes.
- Do not use the device in a humid or damp area, or in water (e.g. bathtub) as it may potentially lead to an electric shock or failure.
- Make sure to keep the device dry during use and storage, as it is vulnerable to moisture. (Please check the humidity level when charging the device.)
- Make sure to use the components provided by the seller
- Do not arbitrarily modify, disassemble or repair the device.

4 PRECAUTIONS

PRECAUTIONS

- You may feel a slight stinging or numbing sensation when using the device in the I mode or L mode. In this case, lower the intensity of the microcurrent to reduce such sensation.
- When using the device in the I mode or L mode, you may temporarily see a flash of light due to the stimulations to the optic nerves.
- If the aforementioned symptom continues even when you are not using the device discontinue using the device.
- Follow the instructions provided in Section 5. Product Care After Use for periodic product care.
- In case of device failure or problems, contact the seller or the local distributor in your country for an inspection or repair. If you do not know the contact information in your country, please contact your local seller for an inspection or repair.

5 WARNINGS / PRODUCT CARE AFTER USE

⚠ CAUTION

Do not use the device in combination with the following types of medical devices, as serious physical impairments may result in case the medical device malfunctions or in the event of an accident:

- Electrical medical devices that are implanted in the body such as an artificial pacemaker;
- Electrical life-support devices such as an artificial cardiopulmonary system;
- Electrical medical equipment with parts that are placed on the body such as an electrocardiogram;
- Medical devices that get attached to the body such as an insulin pump;

※ WARNING: The device should not be used by the following types of people or on the following areas of the body:

- Individuals who have recently undergone cosmetic surgery, or are currently receiving dermatologic care and treatment
- Areas affected with vasculitis, thrombophlebitis, capillary damage, varicose veins, or any other vascular disease, or areas with a skin rash, inflammation, infected skin or swollen skin
- Areas around the eyes (circular muscles within the infraorbitalrim), the temples, the neck (thyroids), the chest or the crotch
- Individuals with an acute disease
- Pregnant women
- Children under the age of 18
- Individuals who exhibit an abnormal response to electrical stimulation
- Individuals with metal allergies
- Individuals with capillary dilatation
- Individuals with a metal implant or any other implants
- Individuals with a history of epilepsy or seizures
- Individuals who have experienced hemorrhaging after an acute trauma or fracture
- Individuals who has prosthetic device for orthodontics
- Individuals who has metal device in their face after surgery

* Product care after use

After using the product, wipe off Faceballs with a wet tissue and a dry towel to prevent foreign matter from remaining on the device.
(Do not wash with water or liquid)

6 PRODUCT SPECIFICATIONS

• LEBODY FACE

Product name [Model name]	LEBODY FACE WHITE: [LEF-100], DARK PURPLE: [LEF-200]
Manufacturer	GTG Wellness #1701~1702, Bundang Suji U-Tower Knowledge Industry Center, 767, Sinsu-ro, Suji-gu, Yongin-si, Gyeonggi-do, Republic of Korea(Headquarters) T. 82-2-3462-5400 Made in Korea
Rating	3.7VDC, 1200mAh
Dimensions	142mm x 51mm x 41mm
Weight	138g

Broadcasting and Communication Equipment
Certification Registration No.

MSIP-REM-sMC-LEF-100

7 WARRANTY

• Warranty Regulations and Conditions for the System

- ① Warranty period: 1 year from the date of purchase (however, for problems resulting from a cause attributable to the customer, the repair service will be provided for a fee)
- ② Scope of the warranty: LEBODY FACE device and accessories
- ③ The warranty will become void and null if the system has been repaired by the customer without the authorization of the manufacturer.
- ④ The manufacturer is only liable for repairing the product within the warranty period.
- ⑤ The warranty is applicable to all of the LEBODY FACE system and accessories that have been used properly.
- ⑥ The warranty does not cover problems resulting from misuse or accidental damage.
- ⑦ Users must adhere to the requirements provided by the manufacturer in relation to the environmental conditions and power rating for use.
- ⑧ The manufacturer cannot be held liable for any problems resulting from arbitrary modification, disassembly and / or repair performed by the end user.

Product Category	Skin Care Device				
Product Name [Model Name]	LEBODY FACE WHITE: [LEF-100], DARK PURPLE: [LEF-200]				
	<table><tr><td>Warranty Period</td><td></td></tr><tr><td colspan="2">※ The internal battery warranty period is 6 months</td></tr></table>	Warranty Period		※ The internal battery warranty period is 6 months	
Warranty Period					
※ The internal battery warranty period is 6 months					
Lot No. /Manufacturing date					
Buyer					
Buyer's Address					
Buyer's Contact Number					

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