




CREATIVE AND INNOVATION ORIENTATED PRODUCTION



Mormed is not just an equipment manufacturer,
is a partner to help you grow and develop your business.

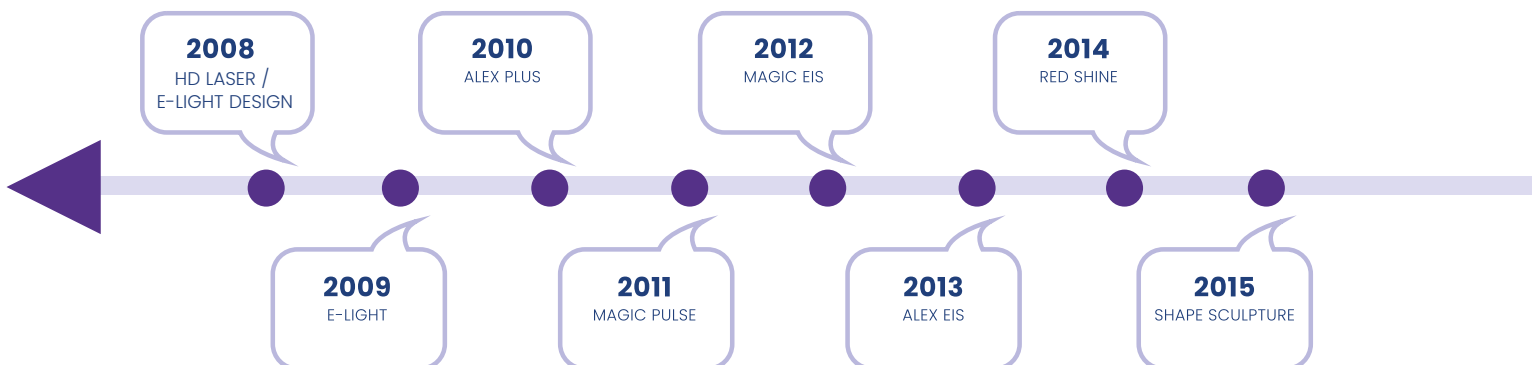

ABOUT US

Our story started in 2004 when Dr Hüseyin Dağar stepped into the field of medical aesthetic medicine. In 2005, it gave its first fruit with the beauty clinic and in 2007, it gave its second fruit with Özlem Dağar's beauty academy. In 2008, we are proud to have led a first in Turkey in the production of equipment by producing a domestic IPL device.

In 2011, while representing Spanish cosmetics and mesotherapy products in Turkey, it added slimming and skin care devices to its range of devices. In this context, it aims to fully meet the need for the equipment and experience required by the salons.

Mormed group has not hesitated to transfer its experience to its stakeholders while continuing its activities in the fields of beauty salon, clinic management, education, production and import.



Mormed's basic understanding is to continue to grow by focusing on winning its customers. It adopts its customers as its own partner. By being with them at every step, it has a service understanding that sees their problems as their own problems and knows their gains as their own gains. This is the only source of advertisement and growth elixir since its establishment.

We aim to make a serious economic contribution to our homeland where we live and produce, and to be among the reputable organisations of our country. By transferring the accumulation and experience of the years we have left behind to the technology we produce, we have continued our work with great seriousness for the development of our country and we continue to continue with all our motivation.

Now we are proud and happy to be the leader of the medical aesthetic sector in our country. Mormed pays its moral debt to this country by selling its products and brands all over the world through its representative offices.


Mormed brand Alex, Magic Eis, Red Shine, HydroOxyPeel, Shape Sculpture, Hobraz, Tatouage, Diode Cardinale, Striaoff, Magic Sound, Big Shaper has gained a place in the globalising world and is moving forward in line with the goals it has set without slowing down. It continues to add added value and employment to the sector with the different technologies it produces by accelerating its R&D studies day by day. We will continue to produce different, innovative, value-adding systems as a result of our work that has been continuing with great devotion for many years.

With the awareness of responsibility given by the leadership in the sector, as Mormed, we attach great importance to customer satisfaction with both our trainings with our professional trainers and our uninterrupted technical support after sales. In addition, we ensure more effective and less consumable application with our business-friendly devices.

We provide local service to all of Turkey with our branches in Konya, Ankara and Istanbul.

We will be happy to see you as a member of this wonderful family.

Dr. Hüseyin Dağar
Özlem Dağar

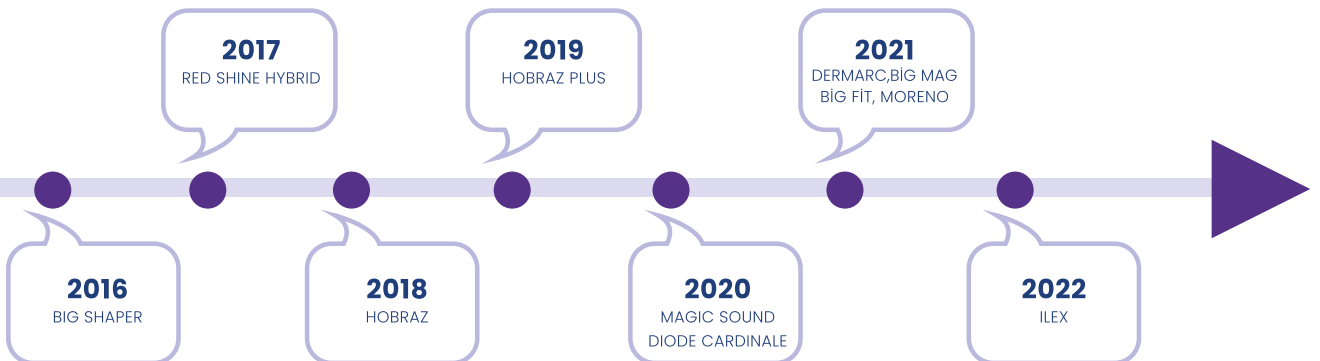

MORMED

Beauty Factory



CONTENTS

Our Company / About Us	02-03
Contents	04
EPILATION DEVICES	
Ilex	06
Diode Cardinale	07
Red Shine Hybrid	08
Red Shine RF+	09
Red Shine	10
Red Shine Fein	11
Alex Eis	12
Alex Premium	13
Alex Plus	14
Magic Pro	15
Magic Eis	16
Magic Pulse	17
Needle Epilation	18
SKIN AESTHETICS	
Tatouage	19
Hydroxypeel	20
Dermarc	21
Strivak - Striaoff	22
Viso Care	23
HD Focus / Eco Hifu / Fraxi HD	24
Aqua Peeling / 4+1 / 8+1 / 12+1 / Skin Care	25
REGIONAL THINNING	
Big Shaper	26
Magic Sound	27
Big Fit	28
Shape Sculpture	29
Hobraz	30
Hobraz Plus	31
Moreno	32
Royal Lift	33
Big Mag	34
Estcultor	35
Bindex	36
Passive Gymnastics/ G5 /Lenf Drenaj	37



NEW GENERATION ALEXANDRITE



Ilex is a laser hair removal device with Alexandrite and Nd:Yag laser technologies with 755+1064 nm wavelength produced by Mormed Group. With its high performance, safe and effective use, it is an ideal choice in hair removal to get the best results expected in terms of both cost and application. Ilex with 755 and 1064 nm wavelengths offers an effective and comfortable hair removal opportunity with the method of focusing on the root of the hair and permanently destroying the hair. Since it creates the cold air itself, there is no extra cooling gas cost.

Principle of operation

755 + 1064 nm wavelength is absorbed by the pigment in the hair root and hair root destruction occurs. Ilex reaches sufficient temperature in the hair follicle in a very short time and allows you to perform a permanent epilation application.

Technology Content

Ilex has laser hair removal technology equipped with a 755 and 1064 nm wavelength head and integrated cold air system. Ilex allows you to perform the procedure in a very short time thanks to its ability to go up to 3 Hz.

Advantages for Business

- With its increased speed up to 3 Hz, it provides the opportunity to apply to more people in a faster and shorter time.
- Practical use with ergonomic head design.
- Saves time thanks to fast and practical application.
- Provides comfortable and continuous sessions thanks to its powerful cooling system.
- Offers ease of use with its ergonomic design and smart digital touch screen.
- It provides an extremely low risk of side effects with a qualified practitioner.
- Success targeted results with effective treatment on all skin types and hair follicles.
- Provides comfortable epilation with low pain threshold.
- It offers 4 seasons epilation with its new generation technology.
- Shot cost is suitable.
- No gas requirement with integrated cold air system.

Benefits to the End User

Even after the first session, the effect is felt epilation application.
Painless and comfortable epilation opportunity.
High level comfort during the whole body epilation process.
Effective epilation in all seasons.
Fast epilation process and short session time.
Possibility to continue daily life after the procedure without interruption.

TECHNICAL SPECIFICATIONS

DIMENSIONS	59 * 120 * 83 CM
WEIGHT	91 KG
ELECTRICAL OPERATION VALUES	220 - 240 V / 50 - 60 HZ
SCREEN	10.4" TOUCH SCREEN
LANGUAGES USED	TURKISH, ENGLISH (OPTIONAL)
NUMBER OF TITLES	1
WORK TECHNOLOGY	ALEXANDRITE LASER - ND: YAG LASER
COLOUR	GREY (OPTIONAL)
HEADINGS USED	ILEX CAP





LIGHTWEIGHT, FAST, PAINLESS
AND NOW, MORE ECONOMICAL
DELIVERING STRONG EFFECT
WITH A LIGHT TOUCH



DIODE CARDINALE 808 is a new generation epilation device developed by Mormed Group, which offers epilation opportunity by delivering the energy generated in the device and this energy to the head through fibre optics instead of the light generated in the classical diode head.

Principle of Operation

CARDINALE 808 is a system that allows the 808nm light inside the device to reach the skin with the help of fibre. While the working temperature was 25-35 degrees in the old systems, with this technology, working comfort is provided with a higher temperature (40-45 degrees). Since the energy is given evenly, the pain sensation is felt at the minimum level, thus increasing the efficiency of the procedure.

Advantages for Business

- Shot cost is more economical than other diode heads.
- Saves time thanks to fast and practical application.
- Less need for technical service thanks to long-lasting and durable consumables.
- Provides comfortable sessions thanks to its advanced cooling system.
- Provides 4 seasons epilation opportunity.
- Ease of use with its ergonomic design, lightweight head and smart digital touch screen.
- It provides an extremely low risk of side effects with a qualified practitioner.

TECHNICAL SPECIFICATIONS

DIMENSIONS	60 * 118 * 73 CM
WEIGHT	60 KG
ELECTRICAL OPERATION VALUES	220 - 240 V / 50 - 60 HZ
SCREEN	10.4" TOUCH SCREEN
LANGUAGES USED	TURKISH, ENGLISH (OPTIONAL)
WAVESIZE	755 - 808 NM (OPTIONAL)
COOLING SYSTEM	CHILLER - AIR - WATER
WORKING SYSTEM	SERIES, SINGLE
COLOUR	GOLD (OPTIONAL)
HEADINGS USED	CARDINALE HEADBOARD

The advantages of the fiber head over the standard diode head

Ultra With its lightweight head (350 gr), it allows epilation without decreasing comfort in prolonged procedures. Since the energy value is higher than other diode technologies, hair destruction is at a high level. Since the energy is given through the device, the head life is longer. Since the energy is given evenly, the pain level is minimised and provides a pleasant and comfortable epilation opportunity. The shot cost is less. The energy reaches the hair root evenly.

Benefits to the End User

Painless and comfortable epilation opportunity.
Safe epilation with a qualified practitioner.
High level of comfort during the whole body epilation process.
The most ideal technology for thick hair on the male back.
Effective epilation in all seasons without postponing your social life.
Fast epilation process short session intervals.



REDShine Hybrid

WITH 3 DIFFERENT WAVELENGTHS

EFFECTIVE, FAST, COMFORTABLE AND NOW
MORE POWERFUL

3 DIFFERENT EXPERIENCES IN ONE
DEVICE



REDSHINE HYBRID DIODE is a diode laser epilation device developed by Mormed Group, which combines the 3 most effective and safest wavelengths (808nm, 755 nm and 1064 nm) for the hair follicle during epilation, allowing simultaneous application.

Principle of Operation

REDSHINE HYBRD technology provides practical and effective hair removal by sending 2 (808nm-1064nm) or 3 (755+808+1064nm) different wavelengths in the same beam to the skin in a single shot and adjusting it according to the hair type. Each of the 3 different laser wavelengths targets a different anatomical structure and different tissue depths in the hair structure.

Technology Content

808 nm Diode: This wavelength has moderate melanin absorption. For this reason, it allows safe epilation on dark coloured skin. **755 nm Alexandrite:** This wavelength is the wavelength that gives the hair its colour and is most absorbed by melanin. It is very effective on light coloured skin and dark coloured hairs. **1064 nm NdYAG:** This wavelength is ideal for very dark skin and tanned skin with low melanin absorption. 1064 nm wavelength is the wavelength that can penetrate deepest into the hair follicle. This helps to achieve successful results in the epilation of deeply located hairs such as the back and groin area.

TECHNICAL SPECIFICATIONS

DIMENSIONS	60 * 120 * 53 CM
WEIGHT	63 KG
ELECTRICAL OPERATION VALUES	220 - 240 V / 50 - 60 HZ
SCREEN	10.4" TOUCH SCREEN
LANGUAGES USED	TURKISH, ENGLISH (OPTIONAL)
WAVESIZE	755 - 808 - 1064 NM (OPTIONAL)
COOLING SYSTEM	CHILLER - WATER
WORKING SYSTEM	SERIES, SINGLE
COLOUR	BLACK (OPTIONAL)
HEADINGS USED	RED SHINE HYBRID HEAD

Advantages for Business

- It helps to eliminate the difference in the number of sessions between hair types as it gives different wavelengths at the same time in a single shot.
- Saves time thanks to fast and practical application.
- Thanks to its long-lasting and durable consumables, the need for technical service is low.
- Ultra cooling system provides comfortable and continuous sessions.
- It offers ease of use with its ergonomic design, lightweight head and smart digital touch screen.
- 4 seasons epilation opportunity with three different wavelength technologies.
- It provides an extremely low risk of side effects with a qualified practitioner.
- Minimises control sessions with serial shot technology.
- Provides a comfortable epilation process with a low pain threshold.
- Shot cost is favourable.

Benefits to the End User

Less pain and comfortable epilation.

It is almost impossible to cause skin damage and stains on dark skin.


High level of comfort during the whole body epilation process.

The most ideal technology for thick hair on the male back.

Effective epilation in 4 seasons.

Fast epilation process short session times.

Ability to continue daily life without interruption after the procedure.

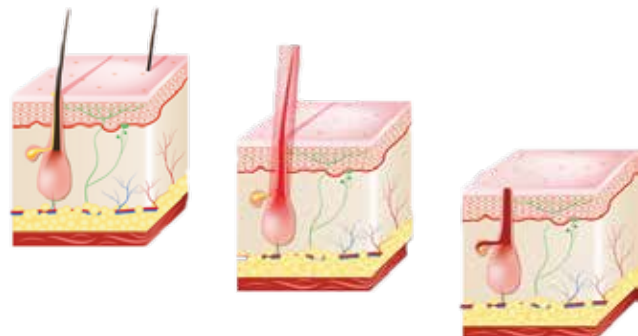

REDShine^{*} +RF

THE MOST EFFECTIVE LASER HAIR REMOVAL SYSTEM FOR FINE AND LIGHT-COLORED HAIRS



REDSHINE DIODE+RF is a new generation radio frequency laser hair removal technology developed by Mormed Group, aiming to rapidly destroy the hair follicles by heating the skin and hair follicles at the highest level thanks to the Diode+RF feature.

Principle of Operation

Thanks to the RF technology, the temperature is first raised up to 45°, enabling more effective epilation with lower energy for the laser. In this way, more efficiency can be achieved on lighter coloured and thinner hair.

Technology Content

RED SHINE DIODE+RF Laser Epilation device with serial pulse system has +RF technology in its head, unlike other serial pulse Diode Laser devices. With radio frequency technology, the degree of heat in the skin is increased during the process and provides much more successful and faster results in epilation. 808 nm Diode, wavelength has medium melanin absorption. Therefore, it enables safe epilation on dark coloured skin. It is in the effective wavelength range in hair follicle destruction.

TECHNICAL SPECIFICATIONS

DIMENSIONS	60 * 120 * 84 CM
WEIGHT	63 KG
ELECTRICAL OPERATION VALUES	220 - 240 V / 50 - 60 HZ
SCREEN	10.4" TOUCH SCREEN
LANGUAGES USED	TURKISH, ENGLISH (OPTIONAL)
WAVESIZE	755 - 808 - 1064 NM (OPTIONAL)
COOLING SYSTEM	CHILLER - WATER - AIR
WORK INTERVAL	0 - 25 J
WORKING SYSTEM	DARBED, SINGLE
COLOUR	DARK PURPLE (OPTIONAL)
HEADINGS USED	RED SHINE -RF TITLE

Advantages for Business

- It is the maximum effective method in fine hair structure and light hair colour.
- With its 808 nm wavelength, it offers safer and faster epilation by reducing the risk of damaging the skin.
- It saves time thanks to its fast and practical application.
- Minimum technical service requirement thanks to long-lasting and durable consumables.
- Comfortable and continuous sessions thanks to the ultra cooling system.
- Ease of use with ergonomic design and smart digital touch screen.
- 4 seasons epilation opportunity with new generation technology.
- Extremely low risk of side effects with a qualified practitioner.
- Provides comfortable epilation with low pain threshold.
- High level hair follicle destruction compared to other diode technologies.
- Less head destruction with extra heating feature.
- More successful and effective results with radio frequency technology.

Benefits to the End User

More effective epilation thanks to Radio Frequency Technology.
Painless and comfortable epilation.
Safe epilation on dark skin.
High level of comfort during the whole body epilation process.
Effective epilation in all seasons without postponing your social life.
Short session times with fast epilation process.
The opportunity to continue your daily life without interruption after the procedure.



REDShine

NOW MORE PRACTICAL WITH 3 INTERCHANGEABLE HEADS ADAPTABLE EPILATION FOR EVERY AREA



RED SHINE 808 DIODE is a fast and effective epilation device with 808 nm wavelength produced by Mormed Group. It provides more comfortable epilation by reducing the risk of damaging the skin with its 808 nm wavelength, which is a longer wavelength compared to other epilation systems.

Principle of operation

REDSHINE technology offers an effective and safe epilation opportunity at the same time by saving time in the epilation process with 10 pulse pulse technique using 808nm wavelength.

Technology Content

808 nm wavelength has moderate melanin absorption. For this reason, it allows safe epilation on dark coloured skin. It is in the most effective and effective wavelength range in hair follicle destruction.

Advantages for Business

- Practical use in narrow and wide areas with 3 different sized interchangeable heads.
- Provides safer and faster epilation by reducing the risk of skin damage with 808 nm wavelength.
- Minimum technical service requirement thanks to long-lasting and durable consumables.
- Saves time thanks to fast and practical application.
- Comfortable and continuous sessions thanks to advanced ultra cooling.
- Ease of use with ergonomic design and smart digital touch screen.
- 4 seasons epilation opportunity with new generation technology
- Extremely low risk of side effects with a qualified practitioner.
- Provides comfortable epilation with a low pain threshold.
- Shot cost is favourable.

Benefits to the End User

Effective and shorter procedures thanks to scaled heads suitable for different skin areas.

Painless and comfortable epilation.

Safe epilation without skin damage and stains on dark skin. High level of comfort during whole body epilation procedures.


Effective epilation process in all seasons.

Short session times with fast epilation process.

Possibility to continue your daily life without interruption after the procedure.

TECHNICAL SPECIFICATIONS

DIMENSIONS	47 * 120 * 47 CM
ELECTRICAL OPERATING VALUES	220 - 240 V / 50 - 60 HZ
SCREEN	8" TOUCH SCREEN
LANGUAGES USED	TURKISH, ENGLISH
WAVELENGTH	755 - 808 - 1064 NM (OPTIONAL)
COOLING SYSTEM	SCREEN COOLED
OPERATING SYSTEM	PULSE, SINGLE
COLOR	GRAY (OPTIONAL)
USED HEADS	RED SHINE HEAD

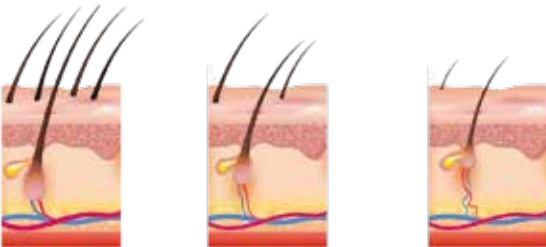

REDShine[®] *fein*

ADVANTAGEOUS COST AND ERGONOMIC DESIGN NEW GENERATION DIODE TECHNOLOGY



REDSHINE FINE Diode is a portable, fast and effective laser epilation device with 808 nm wavelength produced by MORMED GROUP. With its 808 nm wavelength, it offers more comfortable epilation by reducing the risk of skin damage.

Principle of operation

REDSHINE FINE Diode, with its 808nm wavelength technology and 10 pulse pulse technique, saves time in the epilation process and at the same time offers an effective and safe epilation opportunity.

Technology Content

808 nm wavelength has moderate melanin absorption. For this reason, it allows safe epilation on dark coloured skin. It is in the most effective and effective wavelength range in hair follicle destruction.

TECHNICAL SPECIFICATIONS

DIMENSIONS	46 * 38 * 62 CM
ELECTRICAL OPERATION VALUES	220 - 240 V / 50 - 60 HZ
SCREEN	10.4" TOUCH SCREEN
LANGUAGES USED	TURKISH, ENGLISH (OPTIONAL)
WAVESIZE	755 - 808 - 1064 NM (OPTIONAL)
COOLING SYSTEM	WATER - PELTIER
WORKING SYSTEM	DARBED, SINGLE
COLOUR	BLACK (OPTIONAL)
HEADINGS USED	RED SHINE HEADBOARD

Advantages for Business

- Easy to use with its ergonomic and portable design.
- Practical use in narrow and large-scale areas with 3 different sizes of interchangeable heads.
- Safer and faster epilation by reducing the risk of skin damage with 808 nm wavelength.
- Saves time thanks to fast and practical application.
- Minimum technical service requirement thanks to long-lasting and durable consumables.
- Ease of use with ergonomic design and smart digital touch screen.
- 4 seasons epilation opportunity with new generation technology.
- It provides an extremely low risk of side effects with a qualified practitioner.
- Success targeted results with effective treatment on all skin types and hair follicles.
- Provides comfortable epilation with low pain threshold.
- Shot cost is favourable.

Benefits to the End User

Effective and shorter procedures thanks to scaled heads suitable for different skin areas.
Painless and comfortable epilation.
Safe epilation without skin damage and stains on dark skin.
High level of comfort during the whole body epilation process.
Effective epilation in all seasons.
Short session times with fast epilation process.
Possibility to continue your social life without postponing your social life after the procedure.





UNIQUE 3D EPILATION TECHNOLOGY +SHR 2 DIFFERENT EPILATION SYSTEMS COMBINED IN A SINGLE DEVICE



Alex Eis is a double-headed photo-epilation device developed by Mormed Group with its long years of sectoral experience, which offers pulsed radiation technology and Eis cooling head and serial pulse system in a single device.

Principle of Operation

The risk of burning is minimised and the pain sensation is reduced thanks to the pulsed radiation technology and extra cooled head, which includes Alex EIS 610-1200 nm variable wavelength and 3D epilation system. It offers a perfect epilation experience by increasing the speed and comfort with its serial pulse feature.

Technology Content

Epilation can be performed by choosing the ideal system from the 2 different technologies it contains. Alex System: In this system, which can be applied with pulsed radiation technology and square pulse system, effective and safe epilation can be performed even on fine hair and dark skin. In addition, a hygienic process is provided thanks to the personalised silicone head tips that can be used with the head.

+Shr System: While a comfortable application is offered by minimising the feeling of pain with the cooled head, at the same time, much faster and more effective epilation can be performed with the serial pulse system.

TECHNICAL SPECIFICATIONS

DIMENSIONS	50 * 120 * 64 CM
ELECTRICAL OPERATION VALUES	220 - 240 V / 50 - 60 HZ
WEIGHT	70.55 KG
SCREEN	8" TOUCH SCREEN
LANGUAGES USED	TURKISH, ENGLISH (OPTIONAL)
WAVESIZE	610 - 1200 NM
COOLING SYSTEM	PELTIER - WATER - AIR
WORK INTERVAL	10-50 J
WORKING SYSTEM	DARBED, SINGLE
COLOUR	GREY (OPTIONAL)
HEADINGS USED	ALEX HELMET, EIS HELMET

Advantages for Business

- Ease of meeting the technology requests and demands of people thanks to the 2 technologies it hosts.
- Saves time thanks to fast and practical application.
- Ease of use with ergonomic design and smart digital touch screen.
- Provides an extremely low risk of side effects with a qualified practitioner.
- 4 seasons epilation opportunity with new generation technology.
- Shooting cost is appropriate.
- Hygienic application opportunity with personalised replaceable silicone heads
- Effective and comfortable application with 610-1200nm wavelength and 3D photo epilation technology.
- Space and cost saving thanks to its 2 different systems.
- Provides comfortable epilation with low pain threshold.
- Improved cold air system and extra cooled head increases epilation comfort by reducing the risk of skin damage.

Benefits to the End User

Feasible epilation comfort even after the first session with Alex System.

The most ideal technology for those with thick hair structure. Effective epilation in all seasons.

Fast epilation process short session intervals.

After the procedure, you can continue your daily life without interruption.

Painless and Comfortable Epilation opportunity during whole body epilation procedures.

Advantage of using different technology for the same area but different types of body hair.





3D EPILATION TECHNOLOGY CAN BE USED FOR 4 SEASONS AND CAN BE APPLIED TO MANY SKIN TYPES



Alex Premium is a double-headed photo-epilation device developed by Mormed Group equipped with the latest technology, which offers Pulse radiation technology that divides the single shot energy into 2-3 or 4 and serial pulse system with Eis cooling head in a single device.

Principle of Operation

Alex Premium, which contains the most effective 610-1200 nm wavelength and 3D epilation system required for the hair follicle, Pulsed radiation technology and extra cooled head minimise the risk of burning and reduce the feeling of pain. It offers a perfect epilation experience by increasing the speed and comfort with its serial pulse feature.

Technology Content

Epilation can be performed by choosing the ideal system from 2 different technologies.

Alex System: In this system, which can be applied with Pulsed radiation technology and serial pulse system by dividing a single pulse into 2-3 or 4, effective and reliable epilation can be performed even on fine hair and dark skin. In addition, while a personalised hygienic procedure is offered with a silicone head, extra comfort is provided by reducing the pain level with the Air-Liquid-Tec Cold air blowing system.+**Shr System:** While a comfortable application is offered by minimising the feeling of pain with the cooled head, at the same time, much faster and more effective epilation procedures can be performed with the serial pulse system.

TECHNICAL SPECIFICATIONS

DIMENSIONS	60 * 83 * 115 CM
ELECTRICAL OPERATION VALUES	220 - 240 V / 50 - 60 HZ
WEIGHT	90.55 KG
SCREEN	10.4" TOUCH SCREEN
LANGUAGES USED	TURKISH, ENGLISH (OPTIONAL)
WAVESIZE	610 - 1200 NM
COOLING SYSTEM	PELTIER - WATER - AIR
WORK INTERVAL	10-50 J
WORKING SYSTEM	DARBED, SINGLE
COLOUR	DARK PURPLE (OPTIONAL)
HEADINGS USED	ALEX HELMET, EIS HELMET

Advantages for Business

- Saves time thanks to fast and practical application.
- Provides comfortable and continuous sessions thanks to its powerful cooling system.
- Ergonomic design and ease of use with smart digital touch screen.
- Extremely low risk of side effects with a qualified practitioner.
- 4 seasons epilation opportunity with new generation technology.
- Shot cost is appropriate.
- Hygienic application opportunity with personalised replaceable silicone heads.
- Effective and comfortable application with 610-1200nm wavelength and 3D photo epilation technology
- Space and cost saving thanks to its two different systems.
- Improved cold air system and extra cooled head increases epilation comfort by reducing the risk of skin damage.

Benefits to the End User

Ease of technology preference.

Use of a single device for different types of hair on the body.

Painless and comfortable epilation.

Ideal technology for thick hair on the male back.

Effective epilation in all seasons.

Fast epilation process short session times.

High level of comfort during whole body epilation procedures.

Possibility to continue your daily life after the procedure without interruption.

Possibility to epilate without gel.





THE WORLD'S MOST EFFECTIVE 3D AND PULSED EPILATION TECHNOLOGY



Alex PLUS is a new generation photo-epilation device developed by Mormed Group, which provides strong transmission of energy to hair follicles with divisible (2-4) pulse radiation technology, offers successful results even on dark skin and even fine hair, and at the same time allows much more comfortable epilation applications with its integrated cold air system.

Principle of Operation

Alex Plus contains the most effective and effective wavelengths in the range of 610-1200 nm required for hair follicle. Thanks to the pulsed radiation technology and the integrated cold air system, the risk of burning the skin is minimised and the pain sensation is reduced, offering an effective and healthy epilation experience.

Technology Content

Alex System: In this system, which can be applied with pulsed radiation technology, effective and reliable epilation can be performed even on fine hair and dark skin. In addition, while a personalised hygienic procedure is provided with a silicone head, extra comfort is provided by reducing the pain level with the cold air blowing system.

TECHNICAL SPECIFICATIONS


DIMENSIONS	50 * 120 * 50 CM
ELECTRICAL OPERATION VALUES	220 - 240 V / 50 - 60 HZ
WEIGHT	67 KG
SCREEN	8" TOUCH SCREEN
LANGUAGES USED	TURKISH, ENGLISH (OPTIONAL)
WAVESIZE	610 - 1200 NM
COOLING SYSTEM	WATER - AIR
WORK INTERVAL	10-50 J
WORKING SYSTEM	DARBED, SINGLE
COLOUR	WHITE (OPTIONAL)
HEADINGS USED	ALEX CAP

Advantages for Business

- Saves time thanks to fast and practical application.
- Minimises the need for technical service thanks to its long-lasting and durable consumables.
- It offers ease of use with its ergonomic design and smart digital touch screen.
- It provides an extremely low risk of side effects with a qualified practitioner.
- It provides comfortable epilation with a low pain threshold.
- Shot cost is appropriate.
- The ability to process without the use of gel.
- Hygienic application opportunity with personalised replaceable silicone heads.
- Effective and comfortable application with 610-1200nm wavelength and 3D photo epilation technology.
- Integrated cold air system reduces the risk of skin damage and increases epilation comfort.

Benefits to the End User

It provides hygienic epilation with personalised silicone heads.
Painless and comfortable epilation with cold air blowing system.
Safe epilation without skin damage and stains on dark skin.
High level of comfort during the whole body epilation process.
Ideal technology for those with thick hair.
Effective epilation in all seasons.
Fast epilation process short session intervals.
After the procedure, you can continue your daily life without interruption.

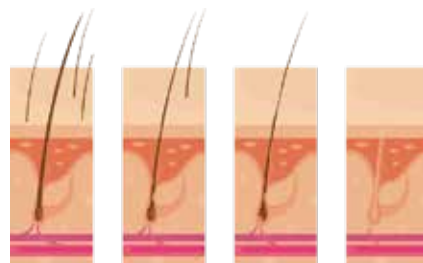



CREATIVE AND INNOVATION BASED
NEW GENERATION IPL TECHNOLOGY
HIGH LEVEL CUSTOMER
SATISFACTION WITH MAGIC PRO



Magic Pro is a photo epilation device with a wavelength of 610-1200nm, which can be shot 10 times per second with serial pulse technology and can be safely processed for 4 seasons with its Eis cooling head. This technology developed by Mormed Group also offers more functional use with its large screen design.

Principle of Operation

Magic Pro, with its rapid-fire technology, can deliver each shot under the skin in milliseconds; It allows you to deliver the right energy to the hair root as much as necessary. Thanks to the extra cooled Eis head, the risk of burning is minimised and the feeling of pain is reduced. With its rapid pulse feature, it offers an ideal epilation opportunity by increasing speed and comfort.

Technology Content

+Shr System: With the extra-cooled Eis head, the pain sensation is minimised and a comfortable application is offered, while at the same time, much faster and more effective epilation can be performed with the Serial pulse system.

TECHNICAL SPECIFICATIONS

DIMENSIONS	60 * 100 * 87 CM
ELECTRICAL OPERATION VALUES	220 - 240 V / 50 - 60 HZ
WEIGHT	55 KG
SCREEN	10.4" TOUCH SCREEN
LANGUAGES USED	TURKISH, ENGLISH (OPTIONAL)
WAVESIZE	610 - 1200 NM
COOLING SYSTEM	PELTIER - WATER - AIR
WORK INTERVAL	10-50 J
WORKING SYSTEM	DARBED, SINGLE
COLOUR	BLACK (OPTIONAL)
HEADINGS USED	EIS HEADER

Advantages for Business

- Saves time thanks to fast and practical application.
- Provides comfortable and continuous sessions thanks to its powerful cooling system.
- Ease of use with ergonomic design and smart digital touch screen.
- Provides hair removal with a qualified practitioner with extremely low risk of side effects.
- Minimises control sessions with +SHR technology.
- Provides comfortable epilation with low pain threshold.
- 4 seasons epilation opportunity with new generation technology.
- Shot cost is suitable.
- Improved extra cooled head increases epilation comfort by reducing the risk of skin damage.

Benefits to the End User

Painless and Comfortable Epilation.
Safe epilation on dark skin.
High level of comfort during whole body epilation procedures.
Effective epilation in all seasons without postponing your social life.
Fast epilation process short session intervals.
After the procedure, you can continue your daily life without interruption.





FAST, SAFE AND EXTRA-COOLED PHOTO-EPIILATION TECHNOLOGY COMFORTABLE EPIILATION



Magic Eis is a photo epilation device with a wavelength of 610-1200nm, which can be shot 10 times per second with +SHR technology and also includes Eis cooling head and high-level serial pulse system. With this technology developed by Mormed Group, it can be safely processed for 4 seasons.

Principle of Operation

Magic Eis can transmit light under the skin within milliseconds with each shot thanks to +SHR technology. Thanks to the extra cooled Eis head, the risk of burning is minimised and the feeling of pain is reduced. With its rapid pulse feature, its speed and comfort are increased, allowing an ideal epilation application.

Technology Content

+Shr System: With the extra cooled Eis head, the pain sensation is minimised and a comfortable application is offered, while at the same time, much faster and more effective epilation can be performed with the serial pulse system.

TECHNICAL SPECIFICATIONS

DIMENSIONS	50 * 120 * 50 CM
ELECTRICAL OPERATION VALUES	220 - 240 V / 50 - 60 HZ
WEIGHT	48 KG
SCREEN	8" TOUCH SCREEN
LANGUAGES USED	TURKISH, ENGLISH (OPTIONAL)
WAVESIZE	610 - 1200 NM
COOLING SYSTEM	SU - PELTIER
WORK INTERVAL	10-50 J
WORKING SYSTEM	DARBED, SINGLE
COLOUR	BLUE (OPTIONAL)
HEADINGS USED	EIS HEADER

Advantages for Business

- Saves time thanks to fast and practical application.
- Provides comfortable and continuous sessions thanks to its powerful cooling system.
- Ease of use with ergonomic design and smart digital touch screen.
- Qualified practitioner provides hair removal with extremely low risk of side effects.
- Provides comfortable epilation with low pain threshold.
- 4 seasons epilation opportunity with new generation technology.
- Shot cost is suitable.
- Improved extra cooled head increases epilation comfort by reducing the risk of skin damage.

Benefits to the End User

Painless and Comfortable Epilation.

High level comfort during whole body epilation procedures. Effective epilation in all seasons without postponing your social life.

Fast epilation process short session times.

After the procedure, you can continue your daily life without interruption.





MAGIC pulse

**MAGIC PULSE NEVER TAKES A BREAK
WITHOUT YOU TAKING A BREAK
FAST, EFFECTIVE AND ECONOMICAL
EPILATION TECHNOLOGY +SHR TITLE**



Magic Pulse is a photo epilation technology that aims to provide a more comfortable and safer epilation application with increased effective process modes. While effective and comfortable application can be performed on many skin types, the application time is minimised by shooting 10 shots per second with serial shooting technology. With this technology developed by Mormed Group, uninterrupted epilation application can be performed for 4 seasons.

Principle of Operation

Magic Pulse, with its Serial Pulse Technology, allows you to deliver each pulse under the skin within milliseconds, while providing sufficient energy to the hair root. With its serial pulse feature, it allows an ideal epilation process with increased speed and comfort.

TECHNICAL SPECIFICATIONS

DIMENSIONS	50 * 120 * 50 CM
ELECTRICAL OPERATION VALUES	220 - 240 V / 50 - 60 HZ
WEIGHT	48 KG
SCREEN	8" TOUCH SCREEN
LANGUAGES USED	TURKISH, ENGLISH (OPTIONAL)
WAVESIZE	610 - 1200 NM
COOLING SYSTEM	PELTIER - WATER
WORK INTERVAL	10-50 J
WORKING SYSTEM	DARBED, SINGLE
COLOUR	PINK (OPTIONAL)
HEADINGS USED	EIS HEADER

Technology Content

+Shr System: While a comfortable application is offered by minimising the pain sensation, at the same time, much faster and more effective epilation procedures can be performed with the serial pulse system. Saves time thanks to fast and practical application. Thanks to its long-lasting and durable consumables, the need for technical service is low. Provides comfortable and continuous sessions thanks to its powerful cooling system. Ease of use with ergonomic design and smart digital touch screen. With a qualified practitioner, the risk of side effects is extremely low epilation process. Provides comfortable epilation for those with low pain threshold. 4 season epilation opportunity with new generation +shr technology. Shooting cost is suitable.

Benefits to the End User

Comfortable epilation with minimum pain
High level comfort during whole body epilation procedures.
Effective epilation in all seasons without postponing your social life.
Fast epilation process short session times.
Possibility to continue your daily life after the procedure without interruption



NEEDLE EPILATION DEVICE

ALL LASER EPILATION AND PHOTO EPILATION

THE MOST EFFECTIVE EPILATION WHEN IT IS INADEQUATE



Needle epilation device is a two-system needle epilation device that includes Flash and Thermolysis systems, which allows application with both classical and painless method. It is one of the permanent hair removal systems that gives the most effective results in permanent hair removal, which is classic but still not outdated.

Needle epilation can be used in most parts of the body, including eyebrows, face, abdomen, thighs, breasts and legs.

Principle of operation

In order to destroy the hair follicle, an electric current is applied into the hair follicle using a very thin needle-shaped electrode or metal probe.

A thin probe follows the hair shaft towards the dermal papillae and at the same time a special electric current destroys the papillae.

The loosened hair shaft is then removed.

Technology Content

Electrolysis, which is a method of removing split hairs from the face or body, is a technology developed to destroy the hair growth centre directly with chemical or heat energy.

Advantages for Business


- Provides effective hair removal in grey, white, red, yellow hair.
- Provides more precise and efficient use with adjustable current pot.
- Provides effective epilation when laser and photo epilation is ineffective.

Benefits to the End User

It offers a more permanent effect rate. It is an effective system in white hair. It increases the success rate in complementary type procedures and Photo Epilation applications.

TECHNICAL SPECIFICATIONS

DIMENSIONS	26 * 34 CM
ELECTRICAL OPERATION VALUES	220 - 240 V / 50 - 60 HZ
POWER	25 WF
WORK TECHNOLOGY	LASH THERMOLYSIS
SCREEN	DIGITAL TIME SETTING 7-SEGMENT DIGITAL DISPLAY
COLOUR	SMOKED
HEADINGS USED	NEEDLE PROBE
OUTPUT FREQUENCY	1.2 MHZ

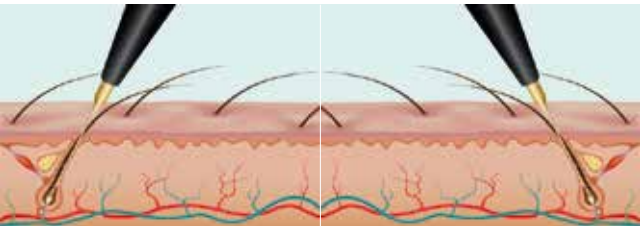

Tatouage

SEGMENTATION & CARBON PEELING QUALITY AND SAFE CARBON PEELING WITH TATOUAGE



Tatouage is a professional Pigment Removal and Carbon Peeling device developed by Mormed Group, equipped with Q-switched & Nd:yag laser technology with powerful energy that covers many colours for successful pigment removal. With 4 different wavelengths, it can be used safely on many skin types. With 532 nm wavelength, you can help reduce red, orange and yellow pigments, with 1064 nm wavelength you can help reduce black and blue pigments. Green coloured pigments can be reduced with 755nm wavelength. Carbon Peeling can be done with 1320nm wavelength. With Tatouage, you can perform eyebrow contour and lip contour beautification procedures, and with the Carbon Peeling head, you can help reduce the appearance of blemishes and shrink the pores.



532nm Wavelength - For light coloured tattoos



755nm Wavelength - For green coloured tattoos



1064nm Wavelength - For dark coloured tattoos



1320nm Wavelength - Carbon Peeling Head

TECHNICAL SPECIFICATIONS

DIMENSIONS	56 * 110 * 54 CM
ELECTRICAL OPERATION VALUES	220 - 240 V / 50 - 60 HZ
WEIGHT	44 KG
SCREEN	10.4" TOUCH SCREEN
LANGUAGES USED	TURKISH, ENGLISH (OPTIONAL)
WAVESIZE	532 - 755 - 1064 - 1320 NM
COLOUR	GREY (OPTIONAL)
HEADINGS USED	TATOUGE CAP

Principle of Operation

With its 4 different heads, it contributes to the correction of skin problems. It produces high quality laser beams with different wavelengths at 532, 755, 1064 and 1320nm, targets pigment particles and the application is made by determining the appropriate wavelength for the pigment colour. Since the beams target only the dye, it can be used safely on many skin types regardless of skin colour. Technology Content

TATOUGE: It is an effective light technology in Depigmentation and Carbon Peeling with its light system that sends energy evenly in pigment removal and head options with a wide range of colours.

Advantages for Business

- Successful and easy depigmentation of various wavelengths Allows you to perform ◦ operation.
- Saves time thanks to fast and practical application.
- Minimum technical service requirement thanks to long-lasting and durable consumables.
- Provides comfortable and continuous sessions thanks to its powerful cooling system.
- Ease of use with its ergonomic design, lightweight head and smart digital touch screen.
- Successful treatments for every colour with 4 different wavelength technologies.
- Extremely low risk of side effects with a qualified applicator.
- Provides a comfortable procedure with a low pain threshold.
- Shooting cost is appropriate.
- Can be used safely on many skin types.

Benefits to the End User

Comfortable Depigmentation process.
Fast and effective process possibility.
Successful Carbon Peeling process with powerful technology.
Advantage of use in all seasons.
Possibility to continue social life without interruption after the procedure.





HydrOxyPeel

INNOVATIVE INNOVATION IN SKIN CARE SYSTEMS ANTI AGING APPLICATIONS



Hydroxy Peel is a new generation skin care device that uses 10-in-one aqua peeling and ultra jet peeling technologies to help you achieve a brighter and more beautiful skin in a comfortable way. The Hydroxy Peel method is a non-invasive system that combines skin cleansing, moisturising and care procedures.

Principle of Operation

Hydroxypeel holistic skin care; It helps to get rid of large or clogged pores, blackheads, uneven skin tone, sebum (oil) imbalance, skin dryness and loss of skin elasticity.

Technology Content

Hydroxy Peel is a multifunctional skin care technology with all the equipment it contains in general skin care.

SOLUTIONS

A - AHA Solution

Glycolic acid (AHA) smoothens the skin by removing dead cells with its peeling feature. At the same time, it provides vibrant and bright skin. It reduces the signs of aging and helps to remove acne and acne scars. It is described as a peeling moisturiser with its water retaining and moisturising effect on the skin.

B - BHA Salicylic Acid Cleansing Solution

Used for oily skin and T-zone. Salicylic acid (BHA) purifies the skin without irritation with its peeling feature. It helps to reduce the problem of oily skin by balancing sebum secretion.

C - Moisturize and Nourish Solution

It is used after procedures with AHA and BHA or on all skin types that need moisture. Thanks to its special formula, it softens, soothes and deeply moisturises the skin.

D - Cleanser

It helps to clean the residues in the system tubes and active system cleaning after each operation.

Advantages for Business

High efficiency with low cost.

Ergonomic design advantage.

Effective applications with heads applicable for each region.

Saving opportunity with the ability to perform many different operations on a single device.

Easy to use thanks to small and lightweight heads.

Benefits to the End User

Ability to continue daily life after the procedure.

Skin care with fast and comfortable procedure.

It offers a quality and effective skin cleansing.

Skin care applications that offer spotless, clean, vibrant, smooth and shine effect with functional use.

TECHNICAL SPECIFICATIONS

DIMENSIONS	41 * 124 * 50 CM
ELECTRICAL OPERATION VALUES	220 - 240 V / 50 - 60 HZ
WEIGHT	55 KG
SCREEN	10.4" TOUCH SCREEN
LANGUAGES USED	TURKISH, ENGLISH (OPTIONAL)
NUMBER OF TITLES	12
WORK TECHNOLOGY	SKIN AESTHETICS
COLOUR	DARK BLUE (OPTIONAL)

LED THERAPY

It is a system that converts light into energy and helps to revitalise the cells without creating a heat effect on the cells by using different coloured light options.

Blue light: Helps to balance the skin's oil production in acne-prone skin.

Green light: Helps to balance the skin tone.

Red light: Helps stimulate collagen synthesis.

Yellow light: It helps to regulate the skin, blood and lymph circulation system and reduce the redness of the skin.

COLD APPLICATION HEAD

Cold application helps to reduce the sensitivity of the skin.

SPRAY

After skin detection, it ensures that the tonic or serum that will meet the needs of the skin is easily spread to the skin.

AQUA PEELING

Aqua Peeling; It helps to achieve a brighter and more beautiful skin in a comfortable way by using water and vacuum methods.

VACUUM

It is used in cases where the skin needs vacuum.

RADIO FREQUENCY

Radiofrequency helps the skin to look more taut thanks to the heat it creates on the skin.

ULTRASONIC PEELING

By creating a series of sound waves on the skin surface, it is a method that helps to cleanse the epidermis from dead cells, blackheads and sebum without irritation.

ULTRA JET PEELING

The air compressor in the device sends air pressure together with water at the time of application. In this way, water and air are sent to the skin.

LIFTING


Lifting helps to train the immobilised facial muscles.

BIPOLAR RADIO FREQUENCY

Radiofrequency helps the skin to look more taut thanks to the heat it creates on the skin.

ULTRASOUND

It helps to ensure the absorption of appropriate and effective products by the skin with sound wave vibrations.

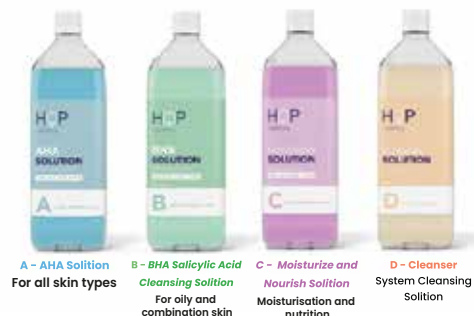

dermarc

CREATE A NEW EFFECT IN
SKINCARE, THE SELF-CARE
OF THE SKIN

THE SYSTEM THAT
ENSURES ITS RENEWAL



What is Dermarc?

The dermarc, which can operate in pulsed and non-pulsed modes, is designed to perform non-contact operation by means of high-frequency current sparks.

In which areas is Dermarc used?

Skin Treatments / Anti Aging Treatments / Stretch Mark Treatments

How does Dermarc regenerate the skin?

It repairs the skin with the thermal effect formed on the skin surface with the current it gives. It triggers collagen and elastin fibres and starts the regeneration process. Dermarc is an effective skin repair system where you can achieve the desired smooth skin on problematic skin tissue.



ARC TOUCH IN SKIN CARE 2 YEAR WARRANTY

TECHNICAL SPECIFICATIONS

DIMENSIONS	27 * 30 * 15 CM
ELECTRICAL OPERATION VALUES	220 - 240 V / 50 - 60 HZ
WEIGHT	4 KG
CURRENT	0.25 A
POWER	PMAX: 55 W
WORKING SYSTEM	ARC SYSTEM
COLOUR	GOLD (OPTIONAL)
HEADINGS USED	DERMARC CAP





INITIATES THE SKIN REGENERATION PROCESS THROUGH VACUUM SKIN CRACK REPAIR SYSTEM



Strivak is a professional skin care device specially developed and produced by Mormed Group, aiming to eliminate skin problems with vacuum technology.

Striaoff is a dynamic-based skin crack rehabilitation system with a success rate of up to 80%, which is performed by combining the new generation strivak device with home use & pro use solutions produced with special formulation and applying the appropriate processing technique.

Principle of Operation

During and after the application, the skin regeneration process is initiated by using home use & pro use solutions, using a powerful vacuum system that can be adjusted according to the structure of the skin contained in the strivak device.

Technology Content

Strivak Device uses powerful vacuum technology that can be adjusted according to the structure of the skin.

Advantages for Business

- Since it appeals to all genders (men and women), you can reach a wider customer base.
- You can make a comfortable application with this skin care system that has almost no complications.
- It is a professional method that gives up to 80% results even in the first session.
- You will be one step ahead of your competitors by increasing your service portfolio.
- You can increase customer satisfaction by providing a solution to this generally sensitive problem.
- Ergonomic and portable design advantage.
- Easy and practical use thanks to its small and lightweight head.
- possibility
- High efficiency with low cost.

Benefits to the End User

Short procedure time with fast and effective application.
Increases skin nutrition with the products used.
Provides psychological and physical relief.
Possibility of procedure without any surgical intervention (non-invasive).
Possibility to continue daily life after the procedure.

TECHNICAL SPECIFICATIONS

DIMENSIONS	41 * 124 * 50 CM
ELECTRICAL OPERATION VALUES	220 - 240 V / 50 - 60 HZ
WEIGHT	10 KG
WORKING SYSTEM	VACUUM
COLOUR	WHITE (OPTIONAL)
NUMBER OF TITLES	1
HEADINGS USED	STRIVAK CAP

SOLUTIONS

PROFESSIONAL USE

NUMBER 1: Hydrogen peroxide is used to disinfect the skin before application to remove germs and bacteria.

NUMBER 2: With the combination of the active ingredients, blood circulation is accelerated and collagen and cell renewal of the dead skin is initiated.

NUMBER 3: With the herbal oils that are beneficial to the body, the effectiveness of the process is increased by using it to slide the head during the process.

NUMBER 4: It is used to accelerate the re-formation process of the skin to remove toxins in our body (poison / harm seen in the organism).

NUMBER 5: It is used to continue the cell regeneration we have restarted, to tighten the skin and to equalise the striae appearance with the skin tone with the beneficial acids in its content.

ZERO: It is an emergency lotion used to neutralise the process and relax the body in case of any unexpected effects.

HOME USE

It is a home care series used to care for stretch marks after a professional stretch mark rehabilitation treatment.

NUMBER 1: By creating a thermal effect on the body, the body is activated by accelerating blood circulation and metabolism and helps to increase collagen production.

NUMBER 2: Moisturises the skin and balances its PH. It helps the skin to breathe and excretion of dead skin from the body.

NUMBER 3: It cleanses the skin from toxins and helps to equalise the skin tone and tighten the skin.

NUMBER 4: It supports collagen production in the skin and helps to maintain cell renewal. Together with these, a smoother and more vibrant skin is targeted.

ZERO: It is an emergency lotion used to neutralise the process and relax the body in case of any unexpected effects.





TECHNOLOGY TO HELP YOU FOR APPROPRIATE TREATMENT SKIN ANALYSER



Viso Care is a skin analyser device used for imaging the subcutaneous skin with the aim of contributing to the diagnosis before care and measuring the success rate of the procedure after treatment with the magic mirror system by using digital image feature thanks to its current technologies such as RGB and UV.

Principle of Operation

Viso Care creates a detailed digital image of the face with the contribution of the magic mirror system and the lighting equipment inside and the simultaneous operation of many high-resolution cameras. In this way, many analyses such as blemishes, wrinkles, skin moisture, sebum and pigmentation are performed.

Technology Content

The skin scanner analyser system is advanced image analysis technology that uses digital image systems through RGB and UV spectrum.
HD Image Quality - ANDROID System - UV Light Technology - RGB Image Quality

TECHNICAL SPECIFICATIONS

ELECTRICAL OPERATION VALUES	220 - 240 V / 50-60 HZ
DIMENSION	42 * 36 * 56 CM
WEIGHT	8 KG
SCREEN SIZE	10.1 INCH
IMAGE QUALITY	HD
COLOURING QUALITY	RGB
WORKING SYSTEM	SKIN ANALYSIS
ANALYSIS SPEED	20 SN
LANGUAGES USED	TURKISH
COLOUR	WHITE
WORK TECHNOLOGY	UV RAY
CORRECTNESS RATIO	%95
JAWS, PEDIMENT MATERIAL	SOFTSILICON

Advantages for Business

- Patient satisfaction and trust can be gained by analysing before diagnosis and after treatment and measuring the success of the procedure.
- More accurate product recommendation can be made with skin analysis.
- Successful results can be presented with pre-application analysis detection.
- Determination of appropriate treatment before and after the procedure.
- Facial images of the client can be compared before and after skin care.
- The success of the treatment can be measured by photographs taken before and after the procedures.
- Ergonomic and portable design advantage.
- High efficiency with low cost.

Benefits to the End User

Fast and comfortable procedure. Increasing the success rate of the procedure to be performed by ensuring the correct diagnosis.



SUCCESSFUL RESULTS WITH PRE-APPLICATION ANALYSIS DETECTION



IDEAL PRODUCT RECOMMENDATION WITH SKIN ANALYSIS



DETERMINATION OF APPROPRIATE CARE BEFORE AND AFTER THE PROCEDURE



WITH ULTRASON TECHNOLOGY TIGHTENING POCKET



ECO HIFU

NO SURGERY FACE AND SKIN LIFT



FraxiHD

MICRO NEEDLE AND RF INCREDIBLE HARMONY OF TECHNOLOGIES



HD Focus 4D is a skin care device that can perform skin tightening with ultrasound waves reaching the subcutaneous tissue without any surgical procedure with its advanced ultrasound technology and different heads for each region. HD Focus 4D offers the advantages of face and skin lifting technique, a popular skin care method, with a comfortable and effective experience without the need for any surgical procedure (non-invasive).

Principle of Operation

With HD Focus 4D New Generation technology, it reaches the subcutaneous tissue and focuses the heat energy directly to the applied area with the ultrasound system. In this way, it stimulates the lower layer of the skin and enables it to regenerate. As a result, it improves the skin texture, reduces sagging of the skin and provides a smooth and aesthetic appearance.

Technology Content

HD Focus 4D has a new generation of 4-dimensional ultrasound technology with different heads with 4 MHz frequency ratio suitable for each region.

3.0 mm head - The whole face except around the lips.
4.5 mm head - Lower cheek and double chin.
6.0 mm head - Breast area and lift.
13 mm head - Under arm, calf and outer calf.
16 mm head - Arm sagging and inner and outer thighs on the legs.

Advantages for Business

◦ Provides high efficiency with low cost. ◦ Ergonomic and portable design advantage. ◦ Provides wide area application with interchangeable heads. ◦ Easy to use thanks to its small and lightweight head. ◦ Creates lifting effect even in a single application. ◦ Effective application advantage with heads applicable for each area.

Benefits to the End User

Skin tightening without any surgical intervention (non-invasive). - Possibility to continue daily life after the procedure. - Fat and cellulite removal with fast and comfortable procedure. - Pore tightening without any surgical intervention. - Helps skin rejuvenation with fast and comfortable procedure. - Improves quality of life with skin tightening feature.

Eco HIFU, odaklanmış ultrason teknolojisi ve her bölge için özel başlıkları ile ameliyatsız yüz ve cilt germe uygulaması yapabilen bir cilt bakım cihazıdır. Cildini gençleştirmeyi, daha genç görünmeyi ve yaşlanmanın cilt üzerindeki etkilerini azaltmayı hedefleyen, yüz germe ameliyatı için kendini hazır hissetmeyen ancak bu tür isteklerine non-invaziv ve girişimsel olmayan bir çözüm arayan kişiler için ideal bir seçimdir.

Principle of Operation

With the technology contained in Eco HIFU, it focuses the heat energy to the subcutaneous tissue with ultrasound waves to the area to be treated and directly stimulates and regenerates the lower layer of the skin. As a result, it improves the skin texture, reduces sagging of the skin and provides a younger and rested skin appearance. Eco HIFU provides skin tension by increasing collagen synthesis. In this way, natural face and skin lifting is provided without recovery time.

Technology Content

Eco HIFU has a new generation ultrasound technology with 5 different heads with a frequency range of 4 Mhz - 7 Mhz suitable for every region. Each cartridge of the Eco Hifu machine reaches more than 10000 pulses.

1.5 mm head - Eye Contour and Forehead
3.0 mm cap - Cheek and Under Chin
4.5 mm head - Under Cheek and Chin
8.0mm head - Arm, Belly and Leg
13 mm head - Arm, Belly and Leg

FRAXI HD, an ideal combination of micro-needling and RF technology, is an ideal device that helps rejuvenate the skin by focusing RF energy precisely directly into a specific skin area with minimally invasive micro-needles. The feature that significantly distinguishes FRAXI HD from laser-based applications is the duration of the procedure and the result. Especially for the face area, it helps to increase skin brightness, as well as reduce acne scars and blemishes. With FRAXI HD, successful results can also be obtained in acne scars, wrinkles, skin blemishes and stretch marks.

Principle of Operation

In the working method of FRAXI HD, it is aimed to open small holes in the skin with micro-needles and help the skin to repair and regenerate them. The radio frequency energy given under the skin heats the skin and with the effect of heating, collagen fibres are shortened and the skin is recovered and tightened. In addition, collagen synthesis is stimulated by the effect of heat and skin repair and regeneration is provided over time.

Technology Content

Fractional RF system is a lifting technology in which micro-needle and radio frequency systems can work simultaneously and stimulate by giving heat under the skin with its heads suitable for different skin depths.

TECHNICAL SPECIFICATIONS

POWER	150 W
CARTRIDES	1.5MM,3.0MM, 4.5MM (FACE) (OPTIONAL)8.0MM, 8.0MM, 10.0MM, 13.0MM., 16.0MM (BODY) (OPTIONAL)
COLOUR	GRAY
CARTRIDGE LIFE	20000 SHOT
FREQUENCY	4MHZ
ENERGY OUTPUT	0.1-2J ADJUSTABLE
VOLTAGE	220/110V;50/60 HZ
PACK SIZE	59*40*25CM
WEIGHT	12 KG

ELECTRICAL OPERATION VALUES	220 - 240 V / 50 - 60 HZ
WEIGHT	31 KG
DIMENSIONS	56 * 36 * 67 C
WORKING SYSTEM	MULTIRASON
COLOUR	GRİ
HEADINGS USED	HD FOCUS BAŞLIK
HEATING TEMPERATURE	65 °C'YE KADAR
PURPOSE OF USE	LIFTİNG
CARTRIDGES USED	1.5MM, 3MM, 4.5MM
WORK TECHNOLOGY	ULTRASOUND

DIMENSION	43 * 35,4 * 26 CM
ELECTRICAL OPERATION VALUES	220 - 240 V / 50-60 HZ
WEIGHT	15 KG
WORKING SYSTEM	SKIN RENEWAL
COLOUR	WHITE
HEADINGS USED	FRAXI HD HEADBOARD
CARTRIDGE TYPES	10/25*64 NEEDLE
SCREEN	8" TOUCH SCREEN
WORK TECHNOLOGY	RADIO FREQUENCY
HEATING TEMPERATURE	UP TO 65 °C



BEAUTY WITH WATER



4+1 SKIN CARE



8+1 SKIN CARE



14+1 SKIN CARE



Aqua Peeling

H2O Aqua Peeling Device is a full-functional skin care device that aims skin care and skin cleansing in general, each equipped with different technology heads (Aqua Peeling, Ultrasonic Peeling, Radiofrequency, Spray, Smooth, Ultrasound Therapy).

Principle of Operation

H2O Aqua Peeling allows pore cleansing, blackhead removal, pore tightening, skin texture and structure strengthening, skin tone equalisation and fine wrinkle removal for an effective and reliable skin care with many different titles in its content.

Technology Content

Aqua Peeling; It is a new generation care application that deeply cleanses your skin and cleanses it from dead skin cells and blackheads, while at the same time providing the moisture your skin needs. The application contains an anti-ageing effect that makes the skin shine radiantly after refreshing care.

Ultrasonic Peeling; It is an effective skin peeling method that allows the skin to be purified from dead cells thanks to a series of sound waves. The resulting sound waves get rid of foreign matter, toxins and cosmetic residues deep in the skin pores.

Radiofrequency; By heating the lower layer of the skin, it helps to tighten and open fine wrinkles, jowl recovery, reduce crow's feet and nasolabial lines, increase collagen production and have a firmer skin.

Spray; Thanks to the tube in its head, it ensures that the products such as tonic, serum, vitamins that we use while skin care are evenly distributed to the skin in a pressurised manner.

Smooth Therapy; Thanks to its powerful cooling engine, it is used to remove the redness on the skin at the end of all the applications you make and to remove the inflammation on the skin by relaxing the skin.

Ultrasound; Ultrasonic sound waves are used to transmit the serums and vitamins that our skin needs to the lower layer of the skin.

Advantages for Business

Increases the variety of procedures in your skin care department with its purchase cost. It provides savings with the ability to perform many different procedures on a single device.

Ergonomic and portable design advantage.

Easy to use thanks to small and lightweight heads.

Benefits to the End User

Ability to continue daily life after the procedure.

Improves the quality of life with its skin tightening feature.

It offers a quality and effective skin cleansing.

Skin care applications that offer spotless, clean, vibrant, smooth and shine effect with its functional use.

TECHNICAL SPECIFICATIONS

ELECTRICAL OPERATION VALUES

POWER

WEIGHT

DIMENSIONS

WORK TECHNOLOGY

COLOUR

SCREEN

HEADINGS USED

220 - 240 V / 50 - 60 HZ

PMAX: 350 W

10 KG

46 * 30 * 43 CM

SKIN AESTHETICS

WHITE

7" TOUCH SCREEN

VACUUM NOZZLE, ULTRASOUND NOZZLE,

SCRUBBER NOZZLE, COLD - HOT NOZZLE,

SPRAY NOZZLE, RF LIFTING

NOZZLEOPSCHOLOGICAL: LED MASK

MULTIFUNCTIONAL SKIN CARE DEVICE

Skin care devices are devices that allow the application of many different techniques required for skin care with their heads with different functions.

Skin Care Device 4+1 Function

1.Steam / 2.Ozone / 3.Loop / 4.Frimator / 5.High Frequency

Skin Care Device 8+1 Function

1.Steam / 2.Ozone / 3.Loop / 4.Frimator / 5.High Frequency / 6.Vacuum / 7.Spray / 8.Iontophoresis / 9.Lifting

Skin Care Device 14+1 Function

1. Ozone / 2. Steam / 3. Illuminated Magnifier / 4. High Frequency / 5. Galvanic / 6. Vacuum / 7. Spray / 8. Spray / 9. Skin Analysis / 10. Frimator / 11. Lifting / 12. LPG Metabolic / 13. Diamond Peeling / 14. Iontophoresis / 15. Hot Hammer / 16. Cold Hammer

Advantages for Business

◦ It is an equipment recommended for use in beauty and aesthetic salons with many functional titles. You can increase the service portfolio and efficiency of your business with this device where you can perform complementary and ancillary applications.

Technical Specifications of the device:

Voltage: 110-220V / Frequency: 50Hz/60Hz / Power: 950 W / Size: 150x 48x 43cm / Weight: 25-35kg

SMALL IN SIZE, BIG IN FUNCTION



BIG SHAPER

REFRESH WITH INNOVATIVE TECHNOLOGY BIG SHAPER GIVES CLEAR RESULTS

Target areas > Arm Area > Abdominal Area
> Cellulite Area > Buttocks Area > Face Area



Big Shaper is a new generation regional slimming device that aims to eliminate cellulite or orange appearance with 4 different technologies such as infrared energy, vacuum, radio frequency and mechanical massage system. It is an ideal choice for shaping, smoothing or toning the body.

Principle of Operation

Infrared energy heats the skin surface and activates the skin's self-renewal mechanism after the heat reaches a depth of 3 mm. Mechanical Massage and Vacuum Technology helps to reduce the bond between the tightened tissues that form the appearance of cellulite by mobilising the subcutaneous tissues. Mechanical Massage process increases the effectiveness of vacuum technology.

Technology Content

Infrared energy and vacuum head supports collagen production by heating the subcutaneous skin and helps to slim down by enabling the fat layer to move with vacuum.

Large head; It helps tightening and slimming by eliminating the appearance of cellulite thanks to its vacuum and mechanical rollers in the leg hip abdominal area.

Small head; It helps recovery and tightening thanks to its vacuum and mechanical rollers in the arm areas.

Face head; It helps to tighten and recover the skin by removing wrinkles with its anti-aging effect.

Advantages for Business

- Combined body care on one platform without taking up much space.
- Painless procedure.
- It offers ease of use with its ergonomic design, lightweight heads and smart digital touch screen.
- With a qualified practitioner, there are almost no side effects.
- It saves time with effective processing in a short time.
- Can be applied in men and women.

Benefits to the End User

Possibility of painless and comfortable application. Ability to continue daily life without interruption after the procedure. Effective application in all seasons. Eliminates the time spent in gyms. Provides rest in the body with massage effect and relieves fatigue.

TECHNICAL SPECIFICATIONS

DIMENSIONS	50 * 122 * 55 CM
ELECTRICAL OPERATION VALUES	220 - 240 V / 50 - 60 HZ
WEIGHT	57.5 KG
LANGUAGES USED	10.4" TOUCH SCREEN
COLOUR	TURKISH, ENGLISH, GERMAN (OPTIONAL)
NUMBER OF TITLES	4
WORK TECHNOLOGY	RADIO FREQUENCY, VACUUM, ULTRASOUND, LED
LIGHT SPECTRUM	850 NM
APPLICABLE REGION	FACE, BODY





2-10 CM IN THE FIRST SESSION WITH MAGIC SOUND OPPORTUNITY TO REVIEW



Magic Sound is a regional slimming device that aims to directly intervene in stubborn regional fats such as Abdominal, Hip, Leg, Arm, Waist and Back by using ultrasonic sound waves, creating a disintegration effect in fat cells with ultrasonic sound waves.

Principle of Operation

Ultrasonic sound waves are sent to the fat tissue with a conductive gel and help to break it down, while at the same time reducing the broken fat cells and providing thinning and tightening in the body, it helps to eliminate the appearance of cellulite to a great extent.

Technology Content

M1 (LARGE) 1.1 MHz Head; It provides thinning and reduction and elimination of fat in unwanted areas.

M2 (SMALL) 40 KHz Head; It provides tightening in thinning and sagging areas without damaging the fat tissue.

TECHNICAL SPECIFICATIONS

DIMENSIONS	41 * 124 * 64 CM
ELECTRICAL OPERATION VALUES	220 - 240 V / 50-60 HZ
WEIGHT	29 KG
SCREEN	10.4" TOUCH SCREEN
LANGUAGES USED	TURKISH, ENGLISH (OPTIONAL)
WORKING SYSTEM	ULTRASOUND SOUND SYSTEM
COLOUR	DARK PURPLE (OPTIONAL)
NUMBER OF TITLES	2
HEADLINE FREQUENCIES	40 KHZ - 1.1 MHZ
HEADINGS USED	40 KHZ HEADER - 1.1 MHZ HEADER

Advantages for Business

- Thanks to the wavelength used, there is no risk of side effects as other tissues are not damaged during application.
- It can be applied in men and women.
- It positively affects the returns by providing a satisfactory result visible from the 15th day.
- Provides permanent results when supported by sports and diet.
- Saves time with short application time.
- It offers ease of use with its ergonomic design, lightweight head and smart digital touch screen.

Benefits to the End User

Painless and comfortable application.

Possibility to continue daily life without interruption after the procedure.

Ease of reaching the desired form in a short time without the need for a surgical procedure.

2-10 cm thinning from the first session.

It offers a permanent solution when supported by sports and diet.



BIGFIT

GOODBYE TO REGIONAL FAT WITH BI FIT



What is Big Fit?

BIG FIT is a non-invasive new generation device developed and produced by Mormed Group and used in the field of aesthetics that targets problematic points that diet and exercise cannot reach by applying ultrasonic sound waves in low and high mode.

Big Fit Working Principle;

The fatty tissue is broken down into small particles that can be easily discharged (excreted). It allows the operator to obtain a constant power with a reduced thermal effect (1-3 MHz and 40 kHz) Therapy (High Frequency Pressure Waves). Thanks to the Pressure Waves that mobilise excess fat deposits, BIG FIT Total Frequency stimulates cellular metabolism by creating an increase in tissue temperature.

TECHNICAL SPECIFICATIONS

DIMENSIONS	41 * 123 * 70 CM
ELECTRICAL OPERATION VALUES	220 - 240 V / 50-60 HZ
WEIGHT	29 KG
SCREEN	10.4" TOUCH SCREEN
LANGUAGES USED	TURKISH, ENGLISH (OPTIONAL)
WORKING SYSTEM	ULTRASOUND SOUND SYSTEM
COLOUR	WHITE (OPTIONAL)
NUMBER OF TITLES	4
HEADLINE FREQUENCIES	:40 KHZ - 1.1 MHZ - 3 MHZ
WORK TECHNOLOGY	ULTRASOUND
APPLICABLE REGIONS	FACE AND BODY

Big Fit Care Before and After;

In order for your applications to be more effective, you should drink plenty of water before and after the application. After the application, a lymphatic massage or a thermal application treatment is recommended to maximise the results.

Application Duration - Applications are 30 minute sessions. You can see results after one session. However, it is recommended in a series of 6-12 sessions.

Pros that distinguish Big Fit from other comparable technologies

It reaches the deepest fat cells with 3 different frequency ranges (40KHZ, 1MHZ and 3MHZ). Thanks to this feature, it provides the opportunity to appeal to 100% of your customer base.

What are the targeted results with the Big Fit application?

Harmonious weight loss
 Reduction in localised fat deposits
 Reduction of cellulite
 Oxygenation of tissues





FOR IMPRESSIVE RESULTS IN
REGIONAL SLIMMING
REGIONAL INSPECTION
TECHNOLOGY PREFERRED
BY AUTHORISED CENTRES



Shape Sculpture is an ideal regional slimming device with its technology that can perform slimming by cooling and dysfunctioning unwanted fat cells in body areas with its specially designed 2 heads, large and small. Basically, it cools the fat in the body and ensures the disruption of fat nutrition.

Principle of Operation

It cools the unwanted excess fat cells formed in different parts of the body and renders them dysfunctional. Fat cells that lose their function are broken down and excreted through lymphatic pathways. This process is completed in a fast, healthy and controlled manner and regional thinning is realised.

Technology Content

The large head non-invasively cools the areas with more fat (hips, belly, buttocks, thighs, calves) and renders them non-functional. The small head non-invasively cools the areas with less fat (arms, back, waist) and renders them non-functional.

Advantages for Business

- It can provide visible effect as it has a reducing effect on fat cells.
- It offers ease of use with its ergonomic design, lightweight head and smart digital touch screen.
- There are almost no side effects with a qualified practitioner.
- It provides ease of application according to the area thanks to its large and small head.
- It is preferred with its ability to be applied in all seasons.
- Can be applied to men and women.

Benefits to the End User

Possibility to continue daily life without interruption after the procedure.

Provides high comfort with painless and painless procedure.

Possibility to reduce fat cells and reach the desired form in a short time without the need for a surgical procedure.

Provides 2-10 cm thinning from the first session.

Possibility of application in all seasons.

Does not damage other tissues during application.

TECHNICAL SPECIFICATIONS

DIMENSIONS	54*120*64 CM
ELECTRICAL OPERATION VALUES	220 - 240 V / 50-60 HZ
WEIGHT	59 KG
SCREEN	10.4" TOUCH SCREEN
LANGUAGES USED	TURKISH, ENGLISH (OPTIONAL)
WORKING SYSTEM	VACUUM - HEATING - COOLING
COLOUR	BLACK (OPTIONAL)
HEADINGS USED	LARGE VACUUM HEAD, SMALL VACUUM HEAD



Hobraz

ACTIVE SOLUTION FOR REGIONAL THINNING RADIO FREQUENCY TECHNOLOGY PREVENTS SAGGING, SLIMS AND TIGHTENS



Hobraz is a regional slimming device developed by Mormed Group by using RF technology in 3 different heads. RF Energy helps to provide collagen production by acting under the skin. In this way, the skin tissue is reshaped and helps thinning and tightening by acting on the fat layer.

Principle of Operation

RF energy helps to increase fat burning and collagen production, while at the same time accelerating blood circulation, giving the skin lost elasticity and tightening it. In radio frequency applications, the water molecules in the cells are activated and reach the fat cells with high water level and help to melt the fat tissue.

Technology Content

Lifting Head; It helps to reduce wrinkles and sagging in the neck, chin, cheek and forehead area. Tripolar Head; It goes deeper into the skin and affects the fat layer and helps thinning.

Multipolar Head; It contributes to tightening by helping to melt the fats in the body and accelerate collagen production.

Business Benefits

Non-invasive application possibility.

It offers a painless and comfortable procedure.

The person returns to daily routine life quickly.

A wide application area with a single device.

It offers ease of use with its ergonomic design, lightweight head and smart digital touch screen.

Minimises the possibility of complications with a qualified practitioner.

Benefits to the End User

Ability to continue daily life without interruption after the procedure.

Provides high comfort with painless procedure.

The desired form is achieved in a short time without the need for a surgical procedure.

Short session times with fast procedure.

RF EFFECT IN SKIN SMOOTHING




TECHNICAL SPECIFICATIONS

DIMENSIONS	60 * 112 * 55 CM
ELECTRICAL OPERATION VALUES	220 - 240 V / 50-60 HZ
WEIGHT	38 KG
SCREEN	10.4" TOUCH SCREEN
LANGUAGES USED	TURKISH, ENGLISH (OPTIONAL)
NUMBER OF TITLES	3
WORK TECHNOLOGY	RADIO FREQUENCY
COLOUR	GREY (OPTIONAL)
HEADINGS USED	LIFTING, TRI-POLAR, MULTITRIPOLAR HEAD

Hobraz plus

LET HOBRAZ PLUS SHAPE YOUR BODY

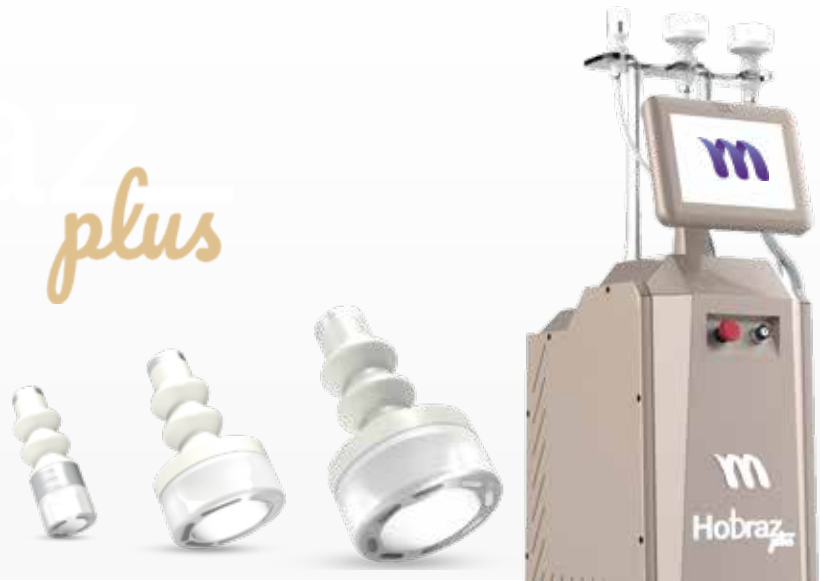

Hobraz Plus, a unique combination of RF energy and Vacuum system technologies, is a new generation regional slimming device that aims to reduce interial fat and smooth the skin.

Principle of operation

RF energy aims to reduce interial fats and accelerates blood circulation with vacuum system technology. In this way, it reduces cellulite and orange appearance, smoothenes the skin and helps the skin to regain lost elasticity.

Technology Content

Hobraz Plus is a technology developed by combining RF energy and Vacuum system, aiming to provide the highest level of process efficiency.

Advantages for Business

- Non-invasive application possibility.
- Offers a painless and comfortable procedure.
- The person returns to daily routine life quickly.
- Wide application area with a single device.
- It offers ease of use with its ergonomic design, lightweight head and smart digital touch screen.
- The possibility of complications is minimised with a qualified practitioner

Benefits to the End User

Possibility to continue daily life without interruption after the procedure.

Provides high comfort with painless procedure.

Possibility to reach the desired form in a short time without the need for a surgical procedure.

It can provide 2-10 cm review from the first session.

Short session times with fast procedure.



TECHNICAL SPECIFICATIONS

DIMENSIONS	53 * 108 * 60 CM
ELECTRICAL OPERATION VALUES	220 - 240 V / 50-60 HZ
WEIGHT	40 KG
SCREEN	10.4" TOUCH SCREEN
LANGUAGES USED	TURKISH, ENGLISH (OPTIONAL)
NUMBER OF TITLES	3
WORK TECHNOLOGY	RADIO FREQUENCY AND VACUUM
COLOUR	GOLD (OPTIONAL)
HEADINGS USED	SMALL, MEDIUM, LARGE HEADER



NEW TECHNOLOGY IN REGIONAL THINNING



What is Moreno?

MORENO is a device used in aesthetic applications that appeals to a new generation wide area (Face and Body) without invasive, produced by Mormed Group using Radio Frequency (RF) energy in Capacitive (CET) and RESISTIVE (RET) 2 modes and VACUUM system.

Deep care with Moreno

Radiofrequency (CET, RET) is much more effective in deep tissues than alternative multipolar radiofrequency systems, safely passing through the body. Besides that, its specific frequency enhances biostimulation and fully respects the natural cell physiology. As a result, the results are not only superficially visible, but also deeply produced with long-lasting effects.

Deep care treatment with Moreno

It gradually and homogeneously increases the temperature in selected areas of the body without disturbing cell physiology, thus providing a comfortable and non-invasive application.

Moreno deep care application session interval

Firstly, good eating habits, adequate fluid intake and light physical activity are recommended for optimum results. A face and neck procedure typically takes 45 minutes, while a décolleté treatment takes about 30 minutes, and a body treatment takes 45 to 59 minutes depending on the area and the individual. 5-6 or 8-10 session intervals may be required depending on the indications and areas of application. The number of sessions is adjusted according to the client's need.

Uses in Moreno aesthetics

- 1- Anti-ageing
- 2- Reduces wrinkles and expression lines.
- 3- Improves the appearance of the eye contour.
- 4- Reduces bags under the eyes.
- 5- Improves the appearance of the double chin.
- 6- Shapes and tightens facial contours.
- 7- Neck and Décolleté
- 8- Body rejuvenation, recovery, reducing the appearance of cellulite.

TECHNICAL SPECIFICATIONS

DIMENSIONS	50 * 120 * 85 CM
ELECTRICAL OPERATION VALUES	220 - 240 V / 50-60 HZ
WEIGHT	44 KG
SCREEN	10.4" TOUCH SCREEN
LANGUAGES USED	TURKISH, ENGLISH (OPTIONAL)
WORKING SYSTEM	RADIO FREQUENCY (CET, RET)
COLOUR	WHITE (OPTIONAL)
NUMBER OF TITLES	4

Moreno deep care result stages

The results of MORENO facial treatments begin to be felt from the first session and become noticeable after the third or fourth application. Body treatments may take longer to show improvement according to BMI (areas where injected fillers are used are not treated). Good eating habits, adequate fluid intake and light physical activity are beneficial for optimum results.

Moreno deep care

Radio frequency is applied with two operating modes: Capacitive mode focuses mainly on superficial soft tissues, suitable for superficial and vascularised tissues. The resistance mode is more suitable for thick, fatty and fibrotic tissues

Are there any restrictions with Moreno deep care?

After the application, you can continue your social activities without having to take any special post-application measures.

By Moreno;

Since the given Radio Frequency energy is harmless to the human body and biological tissues, it has no side effects. However, in some cases, some redness and tingling may be seen on very sensitive skin.

The principle of operation of the Moreno new generation device ;

The face and body treatment uses electrical pulses and heat to stimulate collagen production. The deeper layers of the skin are penetrated by the heat, which can be felt like a deep heat sensation, while the outer layer of the skin is bolstered during the treatment to prevent any damage. New Collagen Radio frequencies are used throughout the body to break down stubborn pockets of subcutaneous fat, which are then gradually eliminated by the body's natural metabolic system. Our MORENO treatment is applied to both the face and body to tighten the skin by stimulating the formation of new collagen. By delivering a controlled heat directly to the tissue that causes collagen contraction and new collagen formation, it helps to reduce fine lines and wrinkles, strengthen the skin and make it look firmer, tighter and healthier. Thus, with this single treatment, your skin receives three important benefits thanks to the combination of different energies and temperatures: Biostimulation (activation of metabolism, collagen and fibroblasts) + Microcirculation (increased metabolism, blood flow and drainage) + hyperactivation (tissue detoxification, fat removal and collagen regeneration). Treatment to remove sagging The treatment can be applied to men and women, both on the face and body, depending on your goal: rejuvenating the face, smoothing wrinkles, revitalising the skin, reducing cellulite, firming the chest, etc. It can also eliminate sagging in the arms and even helps to reduce dark circles and eye bags. The device works with a conductive cream.



NON-INVASIVE BREAST AND BUTT AESTHETICS

6 DIFFERENT HEADS FOR CONVENIENT AND COMFORTABLE APPLICATION



Royal Lift is a new generation skin care device developed by Mormed Group, which aims to solve problems such as sagging and loss of elasticity, allows painless and non-invasive application in the chest and buttocks, and includes Vacuum, Led therapy and Electro therapy (Ems) technologies.

Principle of Operation

Royal Lift; By providing an effective absorption in the breast with vacuum application method, it dilates the blood vessels and accelerates blood circulation. In this way, it increases the rate of breast hormone secretion of the body and helps to add fullness and elasticity to the breast. In problems such as sagging and loosening, it strengthens the chest muscles with the Electrotherapy (Ems) system, helps to tighten and plump the skin by preventing fat tissue with Ems. It provides a fresher and fuller appearance and relaxes with Led Therapy. In addition, for small breasts, it provides a more prominent shape with the 3D pulling force technique by applying a negative pressure.

Technology Content

Royal lift is a skin care technology that aims to lift, plump, tighten and enlarge by using powerful Vacuum, Led Therapy and Electrotherapy technologies with the heads and different adjustment options.

Lifting System: It helps to eliminate deformations such as sagging, loosening, excessive softness that occur in the breasts. Vacuum technology slowly and effectively pulls and stretches the breast tissue so that the stretched cells add elasticity and vitality to the breast tissue.

Plumping and Firming System : Over time, aging, frequent weight gain and loss, childbirth, breastfeeding, improper bra use, etc. It provides a fresh and plump appearance by tightening the breasts that have lost their former appearance for many reasons such as.

Augmentation System: With the scientifically proven vacuum technique for underdeveloped, small and flat breasts, it expands and clarifies the breast volume with fast and effective results.

Advantages for Business

- Successful results and high customer satisfaction with an effective and powerful vacuum system.
- You can achieve high customer conversion in breast-related discontent, which is sensitive and important for women.
- More precise treatment with gradually adjustable suction power.
- Suitable for all breast structures with six different sized vacuum heads.
- It offers the opportunity to process without scar tissue, stretch marks and deformation of the breast.
- More efficient sessions thanks to the increased effect with Mio lifting pads.

Ergonomic and portable design advantage.

- High efficiency with low cost.
- You can make a difference to your competitors by adding this superior technology to your portfolio, which is not available in every business.

Benefits to the End User

It contributes to the quality of life by eliminating breast-related discontent. It is effective in lifting and lifting the buttocks.

It saves time with its fast application feature and effective process results.

Painless and painless comfortable application.

Possibility to continue the daily flow in your life after the procedure.

Gains self-confidence with the aesthetic appearance it provides.

Visible effect after the first session.



TECHNICAL SPECIFICATIONS

DIMENSIONS	46*38*68 CM
ELECTRICAL OPERATION VALUES	220 - 240 V / 50-60 HZ
WEIGHT	24 KG
SCREEN	10.4" TOUCH SCREEN
LANGUAGES USED	TURKISH, ENGLISH (OPTIONAL)
WORKING SYSTEM	VACUUMING - ELECTRIC CURRENTS
COLOUR	BLACK (OPTIONAL)
NUMBER OF TITLES	6
HEADINGS USED	ROYAL LIFT HEADBOARD

BIG MAG

NEW CONCEPT IN FAT BURNING AND SLIMMING

INCREASE IN MUSCLE MASS
NON-INVASIVE BODY
CONTOURING
REDUCTION IN FAT AROUND
THE ABDOMEN
SUITABLE FOR MEN AND
WOMEN
LONG-LASTING RESULTS



HI-EMT/HIFEM

HI-EMT / HIFEM (High Intensity Focused Electromagnetic, High Intensity Focused Electromagnetic energy. Another name is magnetic muscle stimulation.

HI HIFEM

The only way to strengthen core muscle groups is through physical exercise. HIFEM technology can induce almost 50,000 pulses in a single 30-minute session. That's more than 10,000 sit-ups! These pulses stimulate nerves leading to supramaximal muscle contractions (100% muscle contraction), which cannot be achieved by exercise. The muscles are forced to adapt to the resulting muscle thickening. The muscles become stronger, more durable and increase in size.

What's a Big Mag?

It is a new generation body shaping and muscle building device developed and produced by MORMED GROUP using HI-EMT / HIFEM technology. Big mag technology is a new non-invasive body shaping application that only targets muscle as opposed to fat. On average, fat accounts for 25% of total body weight. However, muscle makes up an even larger portion of the human body. Muscle accounts for 42% in men and 36% in women. This technology mimics the action potential of the nerves that cause muscle contractions with a rigour and intensity that you cannot physically realise through voluntary exercise.

Why HI-EMT

- 1-Penetration depth is 8 cm, which covers the entire nerve network and enables contraction of the entire muscle layer.
- 2- Fat apoptosis and the effect of GG: The muscle exercise it performs can never be achieved by physical exercise.
- 3-US studies have shown that the effect it gives is the best.
- 4-Treatment experience and result is good.

BIG MAG Advantages;

- 1- Simultaneous fat burning with growing and strengthening muscles
- 2- Non-invasive application
- 3- No pain sensation
- 4- there is no process of lysisation
- 5- No consumables costs
- 6- Can be applied to different parts of the body (abdomen, buttocks, legs, arms, shoulders, back)
- 7- Saving time with a short 30-minute session
- 8- Process that does not require physical strength

TECHNICAL SPECIFICATIONS

DIMENSIONS	125 * 43 * 76 CM
ELECTRICAL OPERATION VALUES	220 - 240 V / 50-60 HZ
WEIGHT	48 KG
SCREEN	10.4" TOUCH SCREEN
LANGUAGES USED	TURKISH, ENGLISH, GERMAN (OPTIONAL)
COLOUR	WHITE (OPTIONAL)
NUMBER OF TITLES	2
EXIT DENSITY	0-7 TESLA
WORK TECHNOLOGY	HI-EMT / HIFEM

Curious; BIG MAG RESULTS

As long as a good diet is maintained and weight is kept constant, the visual effects of HI-EMT / HIFEM treatment are well maintained. Studies of clients receiving HIFEM treatment have shown that the muscle retains increased size and strength even after one year. This is in contrast to EMS where the effects are temporary.

Target muscles with BIG MAG

Abdominal muscles

Gluteal muscles in the buttocks Thigh and shoulder muscles

Conclusion :HIEMT/HIFEM is very good at strengthening the muscle core. Studies have also shown increased abdominal muscle growth in clients after a full course of treatment. Overall satisfaction levels are very high for the treatment of all muscle groups, especially in the hip area. There are no restrictions associated with the BIGMAG procedure and patients can continue with all their normal daily activities.

With BIGMAG

Research on high-intensity muscle training has shown that a fat-destroying (lipolytic) reaction takes place in the fat lying next to the contracted muscle. HIFEM treatment has been found to affect subcutaneous fat and muscle simultaneously, increasing fat metabolism in a similar way. In scans of clients following HIFEM treatment, reductions in fat and increases in muscle size were found. These changes in fat reduction are long lasting as long as the client maintains their weight with good diet and regular exercise. Studies have shown a reduction in fat and an increase in muscle around the abdomen.



SNOW



ARM



HIP



LEG

igMAG process

HI-EMT/HIFEM is suitable for both men and women. It is excellent for muscle building and fat reduction in those with a body mass index between 19 and 30. Studies have shown that satisfaction levels are higher in clients with a larger body mass index, although the effects are more visible in slimmer clients. A single BIGMEG session results in 50,000 pulse inductions, resulting in 50,000 'supramaximal' muscle contractions. 'Supramaximal' muscle contractions are contractions of the entire muscle as opposed to just 20% to 30% as in a physical exercise session.

The result is hypertrophy (increase in the size of individual muscle fibres) and hyperplasia (increase in the number of muscle fibres as a whole). Studies have shown an increase in muscle mass after 4 sessions. This is seen in the abdominal muscles including the rectus abdominis as well as the lateral obliques. This is known as rectus diastasis and is very common after pregnancy. It is extremely difficult to reduce this gap with physical exercise. Through the associated fat reduction, general improvement around the waist is also achieved. Similar gains are seen in the gluteal muscles of the hips, shoulders and arms and in the quadriceps and hamstrings.

ESTCULTOR

TIGHTEN UP WITH INTENSE SPORTS EFFECT WITHOUT FATIGUE



Device Description

ESTCULTOR is an electro (ems) therapy device developed by Mormed Group, where the highest level of efficiency is aimed with the perfect harmony of micro and galvanic current energies.

Principle of operation,

ESTCULTOR micro and galvanic current energies are used with 5 different head apparatus. It is more effective where there are regional weight and elasticity problems. It provides the removal of cellulite, renewal and repair of the connective tissues (peptide) of weakened fibres. It accelerates collagen and elastin production with the increase of peptides. It renews the skin structure by increasing blood circulation and accelerating cell metabolism with the pressure given during the procedure. Thus, the healing and activation provided within the cell regulates the deterioration in the orange peel appearance, which is the effect of cellulite, and the deformations are visibly reduced as the tissues in the cellulite area regain normal skin texture.

Application Areas Face rejuvenation

Facial renewal application with ESTCULTOR has visible effects that provide facial recovery, increase elasticity and revitalise the skin structure. It helps the skin between 10-15 per cent to gain elasticity and recovery. The recommended session interval is 2 sessions per week for a total of 8 weeks. Although the effect of the procedure is felt as of the first session, visible effects are seen as of 6 and 8 weeks.

Neck and double chin With the decreasing amount of collagen and elastin with the signs of aging, the fat tissue under the skin decreases and causes sagging in the tissue. With ESTCULTOR, it helps to eliminate these problems by creating an intense sports effect on the muscles.

Abdominal recovery and muscles

ESTCULTOR strengthens the abdominal muscles and increases elasticity. The session interval is 2 sessions per week, 6 sessions in total.

Back area Recovery and tightening

ESTCULTOR micro and galvanic current energies and the nerves in the muscle are contracted at the same time and the back recovery process takes place. Intense sports effect occurs in the body. The pressure given during the procedure accelerates blood circulation. Review between 2-12 cm is observed. The recommended session interval is 2 times a week for a total of 8 weeks.

Butt Lift and lift

With ESTCULTOR, micro and galvanic current energies are applied to the desired area at different wavelengths, causing the muscles to contract, creating an intense sports effect, and the blood circulation is increased by the pressure given during the procedure and the buttocks are shaped. The recommended session interval is 2 times a week for a total of 8 weeks.

Leg, calf and inner thigh tightening and recovery

With ESTCULTOR, micro and galvanic current energies are given at different wavelengths and the muscles are contracted, thus creating an intense sports effect. With the pressure given during the process, blood circulation increases and ATP synthesis of the body is stimulated and new collagen is produced. In addition, with the energy generated during this process, the muscles use glycogen and fat in their own tissues. Thus, regional fat tissue is reduced and compression is provided. The recommended session interval is 2 times a week for a total of 8 weeks.

Cellulite and orange peel appearance

ESTCULTOR stimulates the elimination of cellulite and the regeneration and repair of tired fibres and connective tissue (fibroblasts). It accelerates the production of collagen and elastin synthesis. ESTCULTOR helps the tissues in the cellulite area to restructure normally and smooth the area. With the pressure given during the procedure, blood circulation is accelerated, cell metabolism is stimulated and accelerated and the skin structure is renewed.

Advantages for Business

- Wide application area possibility.
- No consumable costs.
- High satisfaction rate
- Appeals to the general customer base.
- Provides comfortable application to customers.

Benefits to the End User

It is painless and painless.
It pleases with visible results.
The skin is also renewed together with the review and tightening.
Without getting tired, your muscles work with intense sports effect.
It also supports your health.

TECHNICAL SPECIFICATIONS

INTRODUCTION	AC110V-220V 50-60 HZ
POWER	350W
OUTPUT FREQUENCY	1hz-200hz
ENERGY	1-100
LANGUAGES USED	TURKISH, ENGLISH (OPTIONAL)
WORK MODEL	7 DIFFERENT MODELS
COLOUR	WHITE (OPTIONAL)



EASY ACCESS TO THE RIGHT TREATMENT WITH BODY ANALYSIS REPORT BODY ANALYSER



Bindex is a body analyser with body composition analysis functions that tests the impedance and composition of elements in all parts of the human body, providing scientific data on weight control, diet, weight loss, muscle development and nutritional balance.

Principle of operation

It applies high-frequency bioelectrical impedance test principle with the 8-point system touching the electrode method, thus helping to prepare an accurate and personalised treatment protocol by measuring weight, protein, fat, water, waist-hip ratio, etc.

Technology Content

- WHR (Waist Hip Ratio)
- Body Fat Analysis
- Bone Weight
- Obesity Diagnosis
- Basic Metabolism
- Lean Weight Indicator
- Nutrition Assessment
- Weight Assessment
- TVS (Body Fluid Percentage)
- PBF (Percent Body Fat)
- Mineral Salt
- Muscle Analysis
- Bioelectric Impedance
- Body Mass Index
- Target Weight Calculation
- Protein Analysis

Advantages for Business

- Body Analysis Report enables correct treatments.
- Ideal treatment recommendation for personalised values.
- Increasing the efficiency of the treatment applied with body values analysis.
- Advantage of working according to the needs of the body by determining the appropriate treatment for pre-procedure.
- Benefits to the End User
- Obtaining accurate information about the body.
- Gaining motivation with the effect of the treatment.
- Correct weight loss.
- A sense of confidence in the accuracy of the treatment to be performed as a result of the analysis.

**Correct diagnosis provides the correct diagnosis and
correct diagnosis provides the correct treatment.**

TECHNICAL SPECIFICATIONS

DIMENSIONS	118 * 50 * 49 CM
ELECTRICAL OPERATION VALUES	220 - 240 V / 50-60 HZ
WEIGHT	28 KG
TEST RANGE	10-200 KG
WORKING SYSTEM	BODY ANALYSIS
COLOUR	WHITE
HEADINGS USED	ELECTRODE HEAD - HANDHELD DEVICE
WORK TECHNOLOGY	HIGH FREQUENCY BIOELECTRIC IMPEDANCE

BODY VALUE ANALYSES



IDEAL TREATMENT RECOMMENDATION BASED ON VALUES



DETERMINATION OF APPROPRIATE TREATMENT FOR PRE-PROCEDURE



25 DIFFERENT TEST RESULTS



Operation Screen Body Analysis Report



ELECTROSTIMULATION PASSIVE GYMNASTICS



Passive Gymnastics is a technology based on the passive exercise method that aims to break down fat and tighten, which enables the device to activate with radio waves by giving low-voltage electricity to the muscles that cannot work in the body with the help of pads. Since passive gymnastics gives the body the effect of exercising and spends energy, calorie burning occurs due to overwork in metabolism. At the same time, the radio waves sent to the body stimulate fat tissues and help the removal of fats from the body.

Principle of Operation

It is connected to the problematic area with the help of pads and low-voltage electricity is supplied and the muscles are activated. Thanks to the presence of more than one pad, it can apply to more than one part of the body at the same time and work the muscles in different areas simultaneously. With different level options, the most suitable setting can be applied according to the person.

Technology Content

It has 8 outputs and 16 pads. It can be operated manually as well as automatically programmed. It is easy to use and provides excellent customer compliance.

Benefits to the End User

It offers the effect of intense exercise in a single session. Improves general circulation in people who do not exercise.

Advantages for Business

Painless, painless and comfortable application. • You can reach much wider customer masses with the variety of services you offer. • More precise applications with adjustable process level setting. • Easy to use with practical installation possibility. • Functional use in upper and lower extremities. • It helps to spend personalised energy with 4 different modes and electric wavelengths. • It helps cellulite and regional slimming applications with Ultrasound and Radiofrequency.

TECHNICAL SPECIFICATIONS

ELECTRICAL OPERATION VALUES	220 - 240 V / 50 - 60 HZ
WORKING SYSTEM	MESSAGE APPLICATION
COLOUR	WHITE (OPTIONAL)
HEADBOARD FEATURES	8 OUTLETS, 16 PED
USABLE AREAS	BEL, ABDOMEN, HIPS
EXTRA FEATURE	CHEST LIFTING

ELECTRICAL OPERATION VALUES	220 - 240 V / 50 - 60 HZ
WORKING SYSTEM	VIBRATION MESSAGE
COLOUR	WHITE (OPTIONAL)
HEADINGS USED	MESSAGE TITLE
WORK TECHNOLOGY	MECHANICALLY ASSISTED VIBRATION MESSAGE
APPLICABLE REGIONS	BELLY, SHOULDER, ARM, LEG

G5 MACHINE BEAUTY WITH MASSAGE



G5 Machine is a professional vibrating massage device used for tightening and cellulite removal with its vibrating system.

Principle of Operation

Mechanically assisted massage therapy is to provide massage application with vibration. It stimulates the muscular and nervous system and at the same time strengthens the circulation of blood circulation and lymphatic system. It softens the fats, beautifies the cellulite tissue, relaxes the muscles and helps to recover.

Technology Content

G5 is a massage technology that can be applied with different titles for all body parts and contains an effective and powerful vibration system.

Advantages for Business

Possibility to combine with another professional Ultra Sound and RF Machine to achieve even greater effect. More precise applications with adjustable intensity level. More functional use with different shape parts such as belly, shoulder, arm and leg. You can increase your service portfolio thanks to its feature that supports physical therapy. You can reach much wider customer masses with the variety of services you offer.

Benefits to the End User

It helps to soften the fat to shape the body.
It helps to separate fat particles and eliminate honey texture.
It helps to effectively relieve muscle tension and fatigue to relax you.
It contributes to relieve muscle pain by eliminating de-toxin and fatigue.
It helps to prevent muscle fatigue and stiffness and also helps to increase the mobility of the joints.
Helps to reduce the lack of muscle exercises to reduce joint oedema.
Helps reduce cellulite by promoting blood circulation.

PRESSOTHERAPY MACHINE LYMPHATIC DRAINAGE Pressure Massage System



Lymph Drainage is a technology that aims to improve skin quality, circumferential thinning, skin and subcutaneous tightening by pressure massage to the body. Lymphatic and venous circulation in the body is accelerated by pressure and helps to drain the fluid accumulated between the veins. It is a method used in the removal of oedema and cellulite caused by various diseases and also increases the effectiveness of regional thinning and body tightening programmes.

Principle of Operation

The person is dressed in lymph equipment and the feet, legs and abdomen are wrapped, then the lymph drainage device is started. The device provides lymph drainage by applying pressure from the fingertips upwards at regular intervals.

Technology Content

The device, which can perform pressure massage on the body, has 10 different programmes. In each programme, the pressure force ranging between 30 and 120 mm/Hg is programmed according to the patient.


Benefits to the End User

It helps to reduce the feeling of oedema, heaviness and pain, especially in the legs.
Supports weight loss, decomposition and fat reduction.
Helps to improve metabolism by promoting blood circulation.
Supports the removal of weight and cellulite.
Contributes to the regulation of lymph flow.
It helps to recover loosened muscle tissue and body sagging.

NOTES

CONTACT

Web www.mormedgroup.com
E-mail info@mormedgroup.com
Phone +90 850 767 0 999

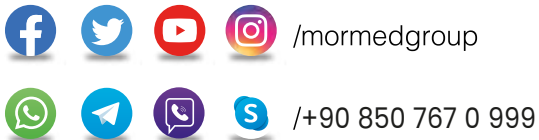

Department Çınardere Mah. Oba Sk. No2/2-3-4-5
Post Code: 34896 Pendik/ İstanbul

Centre Horuzluhan Mah. Kahraman Sk. No 10
Post Code: 42120 Selçuklu/ Konya

Showroom Korkutreis Mah. Sezenler Cad. No: 4/4
Post Code: 06430 Çankaya/Ankara